

# Catalogue



## Contents

|   |            |
|---|------------|
| Shenghua Beauty Catalogue .....   | 1          |
| Company Profile .....   | 1          |
| Cryolipolysis .....   | 错误! 未定义书签。 |
| SH241 Portable Cryolipolysis Slimming Machine .....                             | 3          |
| SH242 Portable Cryolipolysis Cavitation RF Machine .....                        | 5          |
| SH243 Cryolipolysis System with 2 Handles .....                                 | 7          |
| SH244 Super Cryolipolysis with RF Cavitation Device .....                       | 9          |
| SH245 Portable Criolipolisis Equipment .....                                    | 12         |
| SH245B Portable Cryolipolysis Equipment with 2 Handles.....                     | 13         |
| SH245D Updated High Quality Portable Cryolipolysis RF Cavitation Equipment..... | 15         |
| SH246 Cryolipolysis Machine with 3 Cryo Handles.....                            | 17         |
| CRYO6S Portable Cryotherapy Slimming Equipment.....                             | 19         |
| CRV6 Vacuum RF Laser Cryolipolysis Freeze Fat Dissolving Instrument.....        | 21         |
| Fractional RF.....  | 24         |
| NV-V500 Non-Invasive No Needle Fractional RF with 7 Heads.....                  | 23         |
| MR18-25 Cryolipolysis Microneedle Fractional RF.....                            | 27         |
| MR16 Microneedle Fractional RF Beauty Equipment.....                            | 29         |
| Anti-freezing Membrane.....   | 32         |
| SH240 Anti-freezing Membrane/Pad.....   | 32         |
| ETG1 High Quality Anti Freezing Membrane/Pad .....                              | 33         |
| Lipolaser .....   | 35         |
| LS650 Laser Lipolysis Slimming Machine.....                                     | 35         |
| LS651 650nm Diode Laser Lipolysis Machine.....                                  | 36         |
| SH652 Cavitation RF Diode Laser Liposuction Machine.....                        | 38         |
| SH653 Newest 650nm Diode Laser RF Cavitation Equipment.....                     | 41         |
| SH654 Latest Cavitation RF Diode Laser Lipolaser.....                           | 43         |
| SH911 Vacuum Roller Bipolar RF 940nm Infrared Laser.....                        | 44         |
| Vacuum Cavitation RF .....  | 46         |
| SH822 Portable Body RF Face Cavitation Vacuum Roller.....                       | 46         |
| GS6.8 Body Analyse with Vacuum Cavitation Slimming Instrument.....              | 50         |
| GS8.1 Vacuum Cavitation Slimming Beauty Equipment .....                         | 51         |
| GS8.2 E Portable Cavitation Machine .....                                       | 53         |
| M9 Super Body Sculptor.....   | 55         |
| GS9.2 Vacuum Cavitation Slimming Beauty Machine .....                           | 58         |
| GS9.3 Vacuum Cavitation Goethe-shi Electricity Current Beauty Equipment .....   | 60         |
| RU+2 Tripolar RF Cavitation Equipment .....                                     | 61         |
| RU+5 Multipolar RF Vacuum Cavitation Slimming Machine.....                      | 63         |
| RU+6 Cavitation Multipolar RF Tripolar RF Bipolar RF Slimming Machine .....     | 66         |
| SH2.8 Mini Cavitation Vacuum Bipolar RF Multi RF Slimming Machine .....         | 68         |
| SH2.9 Portable Cavitation Vacuum Bipolar RF Tripolar RF Slimming Machine.....   | 71         |

|  |     |
|--|-----|
| Elight IPL RF ND Yag Laser .....   | 72  |
| SH1.1 E-light RF Nd Yag laser Multi-Function Beauty Equipment .....                      | 72  |
| SH1.5 E-light IPL RF Nd Yag Laser.....   | 75  |
| SH1.6 E-light IPL RF Nd Yag Laser Multi-functional Beauty Machine.....                   | 77  |
| SH880 DIY Beauty Equipment.....  | 78  |
| Nd yag Laser .....   | 80  |
| SH3.0 Nd Yag Laser Tattoo Removal Beauty Equipment.....                                  | 80  |
| SH3.1 Portable Q-switched Nd Yag Laser Beauty Machine .....                              | 82  |
| Elight IPL RF .....  | 84  |
| SH4.0 E-light RF Beauty Equipment.....   | 84  |
| SH4.1 Portable E-light IPL RF Beauty Equipment.....                                      | 87  |
| 808 Diode Laser.....   | 89  |
| SH808 808nm Diode Laser Hair Removal Machine.....  | 89  |
| SH810 808nm Diode Laser Hair Removal Equipment.....                                      | 91  |
| SH811 808nm Diode Laser.....   | 93  |
| Asia one 808nm Diode Laser.....  | 94  |
| IPL Device.....  | 98  |
| SH6.0 IPL Permanent Hair Removal Beauty Equipment .....                                  | 98  |
| SH6.1 Portable IPL Beauty Equipment.....   | 100 |
| SH6.2 Portable IPL Photon Machine.....   | 102 |
| SH888 Home Use IPL Device.....   | 103 |
| Infrared Pressotherapy Machine .....   | 108 |
| SH5.0 Infrared Pressotherapy Slimming Machine.....                                       | 108 |
| SH5.1 Air pressure pressotherapy Infrared Equipment.....                                 | 111 |
| SH5.2 Pressotherapy Air Pressure Slimming Equipment .....                                | 113 |
| SH5.4 Portable Air Pressure and Infrared Thermal Equipment .....                         | 114 |
| SH860 Infrared Pressotherapy Slimming Machine.....                                       | 117 |
| RF Radio Frequency .....   | 118 |
| SH7.2 Latest Portable Wrinkle Removal Facial Lifting RF Beauty Machine for Home Use..... | 118 |
| RF3.6 Multipolar RF Vaccum Dissolving Fat System .....                                   | 121 |
| RF395 RF Electric Wave Tighten Skin Instrument.....                                      | 123 |
| Oxygen Injection Instrument .....  | 125 |
| G228A Omnipotent Oxygen Jet Therapy Machine.....   | 125 |
| G882A Hyperbaric Oxygen System.....  | 127 |
| PDT LED Soft Photon Machine .....  | 129 |
| SK8 PDT LED Soft Photon Light Therapy Instrument .....                                   | 129 |
| Fractional co2 Laser.....  | 131 |
| SH9.0 Fractional Co2 Laser Equipment .....   | 131 |
| Breast Enhancer Machine .....  | 132 |
| M7 Super Breast Enhancer machine .....   | 132 |
| N1 Vacuum Breast Enlargement .....   | 135 |
| Micro Needle Roller.....   | 136 |
| BIO50 Micro Needle Roller.....   | 136 |
| Microdermabrasion Machine.....   | 139 |



---

|  |     |
|--|-----|
| NV-60 Diamond Jet Peel Machine.....        | 139 |
| NV1005 Micro Current Facelift Machine..... | 141 |
| Contact US.....                            | 142 |



# Company Profile


## About Us



Shenghua Beauty Equipment Co., Ltd. is an enterprise set up by a group of technical personnel, integrating appearance design, electronic technology development and packing.

Advanced design technique and one-stop service keep customers far from troubles in production and design. Shenghua Company insists on "brand gets upper hand".

From OEM brand-sticking to development and manufacture, we are implementing ISO9001 quality management standard. We promise high quality of all our products, so we have many customers both at home and abroad and get recognized. We help cooperative enterprises set world-class professional image, promote their competitive ability, seize the market and create wealth.

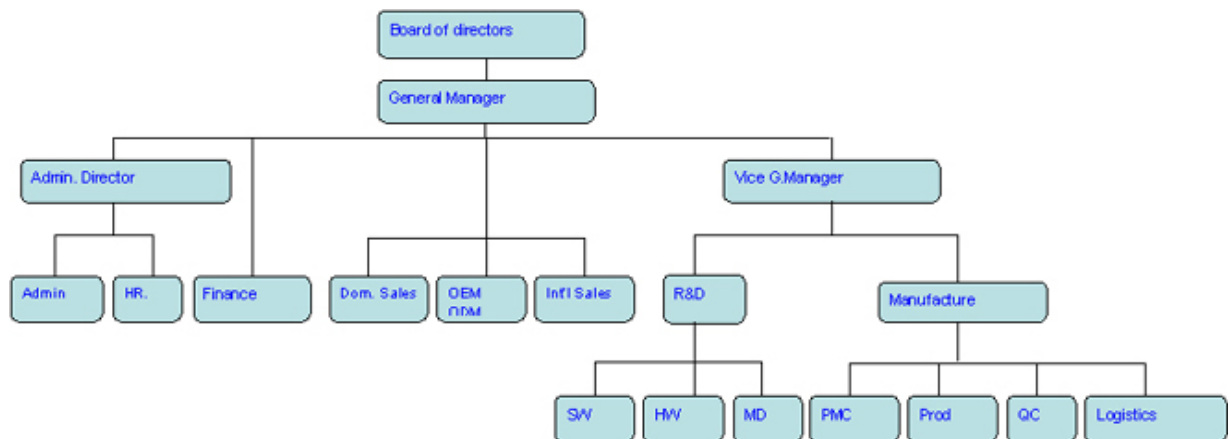

■ 醫療器械生產許可證  
Medical instrument manufacturing license



■ 歐洲CE認證  
European CE certification



■ 產品質量檢測報告  
Product quality inspection report



■ 專利證書  
Patent certificate



# Cryolipolysis

## SH241 Portable Cryolipolysis Slimming Machine



### Specification

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min

Vacuum pressure: 0 - 100Kpa

The host display: 8" TFT color touch screen

The refrigerating temperature: -11 to 5

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Gross Weight: 26kgs

Packing size: 50\*53\*79cm

on both sides of waist), belly and back fat, Cryolipolysis can produce stunning results .

### **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat.

Precisely controlled cooling is applied through the skin to the fat layer, and eliminate local fat. The whole procedure about 1~2 hour.

#### **1. What happens to the fat cells?**

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

#### **2. What will happen after the treatment?**

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

#### **3. How long the effect will last?**

Patients experiencing fat layer reduction show persistent results at least 1 year after the procedure. The fat cells in the treated area are gently eliminated through the body's normal metabolism process. The elimination of fat cells as a result of the procedure is expected to last as long as fat cells removed by invasive procedures such as liposuction.

#### **4. How do the results compare with other procedures such as liposuction?**

liposuction

benefit:

The results of the liposuction procedure are very rapid and dramatic .



Side effect: pain, fatal risk, recovery time even body asymmetry will occur during surgical procedures.

RF thermal therapy or ultrasound

Benefit: can treat deep fat

Side effect: it may hurt the surrounding cells or cause pain

**Cryolipolysis:**

Benefit: Does not harm normal cells

Caution: You must accurately control the position to break down fat and swelling may occur after the treatment.

**Contraindications:**

Liver problems, the Raynaud's disease, severe cold urticaria, cryoglobulinemia, cold fiber hypoproteinemia, severe frostbite, severe diabetes and the elderly, young children, frail and those intolerant of cryotherapy.



## SH242 Portable Cryolipolysis Cavitation RF Machine



### Specification

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min

Vacuum pressure: 0 - 100Kpa

The refrigerating temperature: -11 to 5

Cavitation Power: 150W

Cavitation Frequency: 40Khz

RF Frequency :0.5-10KHz

RF Energy: 50J

The host display: 8" TFT color touch screen

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Packing size: 50\*53\*79cm

### **Body sculpting beauty equipment introduction :**

A number of methods have been attempted for lipolysis including the use of laser, ultrasound, and radiofrequency electric current, since cryolipolysis was developed applying cooling to tissue, was approved as the most effective way to fat reduction. This revolutionary new idea of dissolving fat, which is approved as an important technology progress in fat reduction since it is put into use. For those people who is on diet, do regular exercise but still remove local fat, Cryolipolysis definitely a great gift. For fat intensive parts and less parts, such as love handles (flanks), and back fat. (loose fat above the hip on both sides of waist), belly and back fat, Cryolipolysis can produce stunning results .

#### **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures,

it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat.

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure, there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

#### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.

4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 years warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

## SH243 Cryolipolysis System with 2 Handles



### Parameters

Voltage, AC110V/220V 50-60Hz

Power, < 250W

Temperature , -15°C to 6 °C

Pressure output ,0-100KPa

Packing size: 108.58.69CM

Gross Weight: 55kgs

## **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat. When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat. Precisely controlled cooling is applied through the skin to the fat layer, and eliminate local fat. The whole procedure about 1~2 hour.

### **1. What happens to the fat cells?**

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

### **2. What will happen after the treatment?**

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

### **3. How long the effect will last?**

Patients experiencing fat layer reduction show persistent results at least 1 year after the procedure. The fat cells in the treated area are gently eliminated through the body's normal metabolism process. The elimination of fat cells as a result of the procedure is expected to last as long as fat cells removed by invasive procedures such as liposuction.

### **4. How do the results compare with other procedures such as liposuction?**

liposuction

benefit,

The results of the liposuction procedure are very rapid and dramatic

Side effect,,pain, fatal risk, recovery time even body asymmetry will occur during surgical procedures.

RF thermal therapy or ultrasound

Benefit, can treat deep fat

Side effect, it may hurt the surrounding cells or cause pain

Cryolipolysis

Benefit, Does not harm normal cells

Caution, You must accurately control the position to break down fat and swelling may occur after the treatment

## SH244 Super Cryolipolysis with RF Cavitation Device



### Specifications,

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min

Vacuum pressure: 0 - 100Kpa

The refrigerating temperature: -11 to 5

Cavitation Power: 150W

Cavitation Frequency: 40Khz

RF Frequency :0.5-10KHz

RF Energy: 50J

The host display: 8" TFT color touch screen

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Packing size: 108\*58\*69cm

Gross Weight, 55kgs

Anti freeze membrane, size 28\*32cm

10 pieces , free of charge

### **Body sculpting beauty equipment introduction :**

A number of methods have been attempted for lipolysis including the use of laser, ultrasound, and radiofrequency electric current, since cryolipolysis was developed applying cooling to tissue, was approved as the most effective way to fat reduction. This revolutionary new idea of dissolving fat, which is approved as an important technology progress in fat reduction since it is put into use. For those people who is on diet, do regular exercise but still remove local fat, Cryolipolysis definitely a great gift. For fat intensive parts and less parts, such as love handles (flanks), and back fat. (loose fat above the hip on both sides of waist), belly and back fat, Cryolipolysis can produce stunning results .

### **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures,

it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat.

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure, there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.

2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

## SH245 Portable Cryolipolysis Equipment



### Specifications,

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min



Vacuum pressure: 0 - 100Kpa

The host display: 8" TFT color touch screen

The refrigerating temperature: -11 to 5

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Gross Weight: 26kgs

Packing size: 50\*53\*79cm

### **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat. When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat. Precisely controlled cooling is applied through the skin to the fat layer, and eliminate local fat. The whole procedure about 1~2 hour.

#### **1. What happens to the fat cells?**

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

#### **2. What will happen after the treatment?**

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

#### **3. How long the effect will last?**

Patients experiencing fat layer reduction show persistent results at least 1 year after the procedure. The fat cells in the treated area are gently eliminated through the body's normal metabolism process. The elimination of fat cells as a result of the procedure is expected to last as long as fat cells removed by invasive procedures such as liposuction.

#### **4. How do the results compare with other procedures such as liposuction?**

liposuction

benefit,

The results of the liposuction procedure are very rapid and dramatic

Side effect,,pain, fatal risk, recovery time even body asymmetry will occur during surgical procedures.

RF thermal therapy or ultrasound

Benefit, can treat deep fat

Side effect, it may hurt the surrounding cells or cause pain

Cryolipolysis

Benefit, Does not harm normal cells

Caution, You must accurately control the position to break down fat and swelling may occur after the treatment.

## SH245B Portable Cryolipolysis Equipment with 2 Handles



### Parameters

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min

Vacuum pressure: 0 - 100Kpa

The host display: 8" TFT color touch screen

The refrigerating temperature: -11 to 5

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Gross Weight: 26kgs

Packing size: 50\*53\*79cm

### **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat. When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat. Precisely controlled cooling is applied through the skin to the fat layer, and eliminate local fat. The whole procedure about 1~2 hour.

#### **1. What happens to the fat cells?**

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

#### **2. What will happen after the treatment?**

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

#### **3. How long the effect will last?**

Patients experiencing fat layer reduction show persistent results at least 1 year after the procedure. The fat cells in the treated area are gently eliminated through the body's normal metabolism process. The elimination of fat cells as a result of the procedure is expected to last as long as fat cells removed by invasive procedures such as liposuction.

#### **4. How do the results compare with other procedures such as liposuction?**

liposuction

benefit,

The results of the liposuction procedure are very rapid and dramatic

Side effect,, pain, fatal risk, recovery time even body asymmetry will occur during surgical procedures.

RF thermal therapy or ultrasound

Benefit, can treat deep fat

Side effect, it may hurt the surrounding cells or cause pain

Cryolipolysis

Benefit, Does not harm normal cells

Caution, You must accurately control the position to break down fat and swelling may occur after the treatment.

## SH245D Updated High Quality Portable Cryolipolysis RF Cavitation Equipment



### Specifications,

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min

Vacuum pressure: 0 - 100Kpa

The host display: 8" TFT color touch screen

The refrigerating temperature: -11 to 5

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Packing size: 50\*43\*75cm(L\*W\*H)

### **Beauty equipment introduction:**

A number of methods have been attempted for lipolysis including the use of laser, ultrasound, and radiofrequency electric current, since cryolipolysis was developed applying cooling to tissue, was approved as the most effective way to fat reduction. This revolutionary new idea of dissolving fat, which is approved as an important technology progress in fat reduction since it is put into use. For those people who is on diet, do regular exercise but still remove local fat, Cryolipolysis definitely a great gift. For fat intensive parts and less parts, such as love handles (flanks), and back fat. (loose fat above the hip on both sides of waist), belly and back fat, Cryolipolysis can produce stunning results.

### **Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures,

it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat.

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure, there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

### **Feature**

1. The vacuum Roller including the RF which is much more effective
2. The cavitation head can emit RF at the same time.
3. Two different size of tripolar tips for face and body

4. Cooling probe for after treatment care
5. Big touch screen with easy operation

## SH246 Cryolipolysis Machine with 3 Cryo Handles



### Specifications:

Power: 500W

Vacuum: 650mmHg

Pump flow: 60L / min

Vacuum pressure: 0 - 100Kpa

The host display: 8" TFT color touch screen

The refrigerating temperature: -11 to 5

Working voltage: AC220, 230V 50Hz, AC110 - 120V 60Hz;

Packing size: 50\*43\*75cm(L\*W\*H)

**Beauty equipment introduction :**

A number of methods have been attempted for lipolysis including the use of laser, ultrasound, and radiofrequency electric current, since cryolipolysis was developed applying cooling to tissue, was approved as the most effective way to fat reduction. This revolutionary new idea of dissolving fat, which is approved as an important technology progress in fat reduction since it is put into use. For those people who is on diet, do regular exercise but still remove local fat, Cryolipolysis definitely a great gift. For fat intensive parts and less parts, such as love handles (flanks), and back fat. (loose fat above the hip on both sides of waist), belly and back fat, Cryolipolysis can produce stunning results.

**Cool Technology**

As triglyceride in fats will be converted into solid in particular low temperatures,

it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat.

During the procedure, exposure to cooling causes fat cell apoptosis, which leads to the release of fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process, much like fat from food.

The procedure is non-invasive. Most patients typically resume their routine activities, including returning to work and exercise regimens on the same day as the procedure. Transient local redness, bruising and numbness of the skin are common side effects of the treatment and expected to subside in couple of hours. Typically sensory deficits will subside within 1~8 weeks.

With this non-invasive procedure, there is no need for anesthesia or pain medications, and no recovery time. The procedure is comfortable for most patients can read, work on their laptop computer, listen to music or just relax.

**Feature**

1. The vacuum Roller including the RF which is much more effective
2. The cavitation head can emit RF at the same time.
3. Two different size of tripolar tips for face and body
4. Cooling probe for after treatment care
5. Big touch screen with easy operation

## CRYO6S Portable Cryotherapy Slimming Equipment



### Technical Parameter:

- 1.Voltage: AC110V/220V 50-60Hz
- 2.Power: <500W
- 3.Temperature: -15 to 5
- 4.Output pressure:0-100KPa
- 5.Vacuum: 10-80KPa
- 6.Cool Liquid: Distilled water or special coolant
- 7.Fuse: T3. a5AL250V
- 8.Ambient temperature:5-40
- 9.Relative Humidity: <80%
- 10.Atmospheric pressure: 86Kpa-106Kpa
- 11.Packing Size: 50\*52\*80cm



12. Gross Weight: 28kgs

### **Cryolipolysis Machine Introduction**

A number of methods have been attempted for lipolysis including the use of laser, ultrasound, and radio frequency electric current, since cryolipolysis was developed applying cooling to tissue, was approved as the most effective way to fat reduction. This revolutionary new idea of dissolving fat, which is approved as an important technology progress in fat reduction since it is put into use. For those people who is on diet, do regular exercise but still remove local fat, Cryolipolysis definitely a great gift. For fat intensive parts and less parts, such as love handles (flanks), and back fat. (loose fat above the hip on both sides of waist), belly and back fat, Cryolipolysis can produce stunning results .

### **Cryolipolysis Machine Technology**

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat. Precisely controlled cooling is applied through the skin to the fat layer, and eliminate local fat. The whole procedure about 1~2 hour.

### **Functions:**

Use the low temperature to freeze the fat cell, then to make the fat reduction.

It is the real cryolipolysis/freeze fat technology to make the local part fat reduction.

The effective result can be seen after one time treatment.

This cryolipolysis technology is the best way to reduce the fat ,not reduce the water.

You can see the fat reduce about 26% after one time treatment.

## CRV6 Vacuum RF Laser Cryolipolysis Freeze Fat Dissolving Instrument



### Technical Parameters

1. Input : 220V/ 50Hz 110V/60Hz
2. Power : < 500 VA
3. Heat : 37°C~45°C
4. Cryo: 5°C~-10°C
5. Vacuum: 10~80Kpa
6. Light : Red( 630nm) Green( 570nm) 50mWX
7. Cooling liquid: Pure water or special cool
8. Volume: 44cmX48.5cmX98cm
9. Fuse: T3. 15AL250V
10. Ambient temperature: 5°C~40°C
11. Relative Humidity : ≤80%

12. Out put RF frequency: 0.8MHz to 2.6MHz

### **1.Cryolipolysis Theory**

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat. When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And The fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat. Precisely controlled cooling is applied through the skin to the fat layer, and eliminate local fat.

### **2. Vacuum Multipolar RF Theory**

( Vacuum Multipolar RF +Red Light ) deploys heat using multiple Radio Frequencies in 3 different treatment modes that target different tissue depths, from superficial layers of the skin to the hypodermis. The 4th mode - multi-channel mode - creates a stronger, more accurately distributed thermal response that travels through all skin layers. this new and innovative technology also facilitates easier penetration of RF energy which travels more easily through pre-heated tissue.

### **3.Blue/ Red Light Theory**

Red Light main function: The red light can accelerate the growth of collagen fiber by photochemical biology activity. The red light act on chondriosome in the hypoderm to generate more power which speed up the synthesis of DNA, RNA, and a lot of collagen and fiber is produced and accelerate the elimination of waste or dead cell, then to alleviate the symptom of aging skin, large pore, fine wrinkle, and muddy complexion.

Blue Light main function: The visible blue light emitted by LED can be perfectly absorbed by acne bacillus, and the metabolism of bacillus will generate a mass of singlet oxygen which is harmful to the acne bacillus so as to kill them and eliminate acne in the skin. Photochemical reaction also accelerates the collagen fiber generation.

### **4. Cold Laser Theory**

Cold laser is a pain-free, non-invasive device used for reducing fat by shining the emitted light (laser) over the skin. Cold laser is a revolutionary new fat reduction technology.

A smart alternative to surgery, Cold laser uses a low-level laser to target the stubborn spots that don't respond to regular diet and exercise.

Cold laser gives fast results, with improvements showing immediately after the first treatment and over the course of treatments. Cold laser is ideal for someone who has made a decision to improve their overall diet and exercise regime, or has struggled to get results from a particular area despite regular exercise and healthy eating.

Cold laser is a spot fat reduction technique that help lose weight in the areas that are traditionally much harder to target.

# Fractional RF

## NV-V500 Non-Invasive No Needle Fractional RF with 7 Heads



**Specification:**

| Functions                 | Specifications                    |
|---------------------------|-----------------------------------|
| 1 Fractional RF Frequency | Biopolar & Multipolar 2MHz / 4MHz |

|    |                                   |  |
|----|-----------------------------------|--|
| 2  | Fractional RF <b>Energy Level</b> | 25 Joules  |
| 3  | Fractional RF <b>Handpiece</b>    | 2 PCS  |
| 4  | Fractional RF Tips <b>Types</b>   | 5 Types  |
| 5  | Fractional RF Tips <b>Life</b>    | 3000 Shots   |
| 6  | Fractional RF Tips <b>Size</b>    | 0.38cm <sup>2</sup> , 2.0cm <sup>2</sup> , 3.61cm <sup>2</sup> |
| 7  | Fractional RF Tips <b>Pins</b>    | 24 Pins, 146 Pins*2, 64 Pins, 196 Pins                         |
| 8  | Digital System                    | Digital Auto counting  |
| 9  | Thermal RF Types                  | 2 MHz / 4 MHz  |
| 10 | Thermal RF Frequency              | Bipolar & Multipolar 2MHz / 4MHz                               |
| 11 | Thermal RF Mode                   | Pulse & Continuous   |
| 12 | Safe Protection                   | Coded Lock   |
| 13 | Display                           | 7 inch LCD Touch Screen  |
| 14 | Power Supply                      | AC 100V-120V/200V-240V/50Hz/60Hz                               |

**Treatment programs:**

- 1.Facial Mesotherapy: Sculptures face and chin curve contour, to reconcile skin texture and feel body.
- 2.Sculpting:Firming sagging loose skin of the upper arm; improve stretch
- 3.Marks eye skin rejuvenation: Slow down the sagging of the eyelids, fade away crow's feet and eye fine lines.

**Working principle:**

Fractional RF system offers advanced fractional multipolar radio frequency technology.The treatment remains ablations zones of micron depth including inside deep necrosis regions, it promotes more effective collagen remodeling compare with traditional resurfacing treatment.

**Applications:**

- 1-Wrinkle
- 2-Scar
- 3-Pore
- 4-Hypertrophic texture
- 5-Pigmentation

6-Tightening

7-lifting

8-Acne/Acne scar

9-Contouring

**Technology Advantages:**

1) Easy Operation

Intelligent Color Screen System installed with Voice Prompt System, the operator can use it easily and exactly.

Nenergy of RF plays a role in dermis directly, therefore, there is no heat summation on the dermis to avoid the possibility of blister and pigment precipitation problems.

2) No Pigment

The energy of RF plays a role in dermis directly, therefore, there is no heat summation on the dermis to avoid the possibility of blister and pigment precipitation problems.

3) Less Time

It just takes a few minutes (5—10 minutes) to do the treatment, such as Wrinkle Removal, Whitening,

Shrinking Pores, Acne Scars Removal, Skin Rejuvenation

4) No Side Effece

The recovery time is short. It won't affect the daily life after treatment. You can wash your face and make-up as usually.

5 Two Treatment Ways

There are two treatment ways with Dual Matrix and RF Microneedle to meet the needs of different clients's psychology and physiology.

Compared with the 1550 Fractional Laser in the market, the effect of Floracel is more wonderful, the treatment is safer.

It is warm welcomed around home and abroad.

## MR18-25 Cryolipolysis Microneedle Fractional RF



### Technical parameters:

Power: <80W

Output frequency: 1.6MHz

Output wavelength: PDT 650nm

Freezing temperature: -5°C

Ambient temperature range of: -40°C— +55°C

Relative humidity range: ≤93%

Atmospheric pressure range of: 50kPa-106kPa

voltage: 220V/50Hz, 110VA/60Hz

Package Size: 56\*66\*145cm

Gross Weight: 30 kgs

### Working principle:

Fractional RF needle system through the special design of many fractional needle point array, the high speed digital motor controls order through the epidermis and the dermis accurate to control the depth of 0.25-3 mm, again by the end of lattice needle release RF, stimulate collagen and elastic tissue, and is not thermal damage skin layer potential. When skin tissues through this process, the epidermal layer is safe, RF energy can well to penetrate to the dermis, stimulate collagen protein hyperplasia, not only is the best way to improve the scars, but also long-term tightening skin wrinkles induced the good method. Cryo helps SRF and MRF play its function.

**Feature:**

1. Easy to control the needle into the depth, agile adjustable
2. Superficial fractional RF and microneedle fractional RF hand stick hand bar
3. The allocation of 25 pin, 49 pin and 81 pin replaceable type RF needle head
4. The allocation of replaceable type RF point array work head
5. RF energy fine adjustable
6. MRF+SRF+PDT+CRYO
7. The mature stable eight inches of true color liquid crystal display output and input system
8. Configuration sellers market plan management system

**Function:**

1. Improve pores, scars, deep wrinkles, acne, pigment calm.
2. It makes the skin smooth, tighten pores
3. It has a minimum of pain and provide accurate energy
4. The minimally invasive make even less epidermal damage
5. It makes the RF direct penetrate to the target area, in order to obtain a sufficient effect.
6. Scar removal including acne scar reduction
7. Applies to: stretch marks, anti-aging, anti-wrinkle treatment, hair loss/hair restoration, Hyperpigmentation treatment.

**What is the Liftron?**

48 microneedles can be inserted from the epidermis to the dermis and generate RF energy to stimulate the fibroblasts in the dermis effectively to improve the elasticity of the skin and wrinkles.

**Handpiece**

- Convenience for user
- Minimize pain for treatment
- Convenient control of the depth of the microneedle
- Insulated needles for epidermis without any thermal damage
- Convenient replacement of the Microneedle tip



car goes away quickly

**Effect**

Microneedle(Stamping Effect)

This is perceived by the skin as damage, which triggers the release of growth factors that trigger the production of collagen and elastin.

**RF Effect**

Radio Frequency can be penetrated to dermis without any damage to the epidermis therefore elastic skin can be made in a short time.

**Combination Therapy**

Can be combined fractional RF with another laser treatment at the same day.

**MR16 Microneedle Fractional RF Beauty Equipment**



### **Specifications**

Power Input: 110/220Vac 60Hz  
Dimensions: 47(W) 46(D) 110(H) cm  
Weight: 25kg  
Matrix Size : 0.8x0.8(cm)  
Disposable Tip Matrix : 25 electrode needles  
Output Frequency: 1MHz  
RF Output Energy Power : 25W

### **Effect**

Microneedle(Stamping Effect)

This is perceived by the skin as damage, which triggers the release of growth factors that trigger the production of collagen and elastin.

### **RF Effect**

Radio Frequency can be penetrated to dermis without any damage to the epidermis therefore elastic skin can be made in a short time.

Combination Therapy

Can be combined fractional RF with another laser treatment at the same day.

### **Characteristics**

Easy control of the depth of the needle  
Easy Control of RF output  
No taking people away from work or normal life cycle.  
Short downtime  
No pigmentation, no burning.

### **Treatment Principles.**

It can treat pores, scar, deep wrinkles, acne, pigmentation at the same time.  
It produces collagen and elastin in the skin and improves the wrinkles.  
It makes skin smooth and tightens the pores by RF.  
It has the minimal pain and deliver the energy completely.  
It goes to the target without any damage on the epidermis.  
It makes RF penetrated to the targeting area directly to get the enough effect.

### **Indications**

Scar Removal – including acne scar reduction  
Stretch Marks Removal  
Anti Ageing Therapy  
Anti Wrinkle Therapy  
Hair Loss Therapy / Hair Restoration  
Hyper Pigmentation Treatment

System composition

Main unit

Handprobe(1pc)

Disposable tip(1pc)

Foot Switch(1pc)



# Anti-freezing Membrane

## SH240 Anti-freezing Membrane/Pad



### Specifications:

Big size 28\*32cm

Weight:70g

Protect skin for cryo treatment

Validity 3 years

### Delivery Time:

1-2 days

## ETG1 High Quality Anti Freezing Membrane/Pad



### Specifications:

Big size 28\*28cm

Weight 60g

Protect skin for cryo treatment

Validity 3 years

### Delivery Time:

1-2 days



**Delivery Way,**

Door to door, by EMS, DHL.

To airport, by Aircargo,

To Sea port, by Sea



# Lipolaser

## LS650 Laser Lipolysis Slimming Machine



### 1. Specification

|                   |                         |
|-------------------|-------------------------|
| input voltage     | AC110V/220V 50-60Hz     |
| power consumption | <500VA                  |
| screen            | wide color touch screen |
| Output frequency  | 1-1000Hz 40K 5M         |
| Energy output     | 18x40mW(total 0.72w)    |
| Wave length       | 650nm                   |
| OEM service       | available               |

|                  |  |
|------------------|--|
| Application      | Beauty salon, health center, school, person etc. |
| Features         | Easy to operate, study                           |
| Training service | We provide customers training service            |
| carton size      | 56x53.5x124.5cm                                  |
| Working theory   | Laser, cavitation, rf, vacuum                    |
| N/W              | 20KG   |
| warranty time    | 1 year   |
| delivery time    | 2-3 days   |

## 2. Operation

Low levels of laser energy emitted by LS650 stimulate the fat cell membranes, changing their permeability. The cells lose their round shape, and intracellular fat is released. Then the fatty triglycerides flow out of the disrupted cell membranes and into the functions with no harmful physiological effects, to be used as an energy source for peripheral nerves. It is not merely the liquefaction of fat but instead it is the instant breakdown of the fat cells, otherwise known as lipolysis.

LS650 the intelligent alternative to lipolysis.

LS650 is the non-invasive procedure to effectively remove excess fat.

## 3. Function

LASER, CAVITATION and VACUUM:

Innovative Treatment for reducing Fat and Reshaping the Body .

LS650 provides you a way to succeed in weight loss and body contouring without pain, needles and down time.

## 4. Company advantages:

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.



## LS651 650nm Diode Laser Lipolysis Machine



### Theory:

laser for fat reduction,use the 650nm soft laser to destroy the deep fat cells, broke the cell membrane,change the fat to triglyceride, then leave out the body thanks to metabolism.

The laser heat can Stimulate the subcutaneous dermal tissue collagen and elastic fibers hyperplasia,make the treated area skin tighten and elasticity .

### Parameters:

Output frequency: 1-1000HZ,40KHZ,5MH

Technology: laser

Energy output: 18\*40MW(Total 0.72W)

Wavelength: 650nm

Power: 500VA

Hand piece: x2 laser namd, X4 laser hand

N.W/G.W: 8.5kgs/15kgs

Size: 50cm\*54cm\*33cm

**Treatment procedure:**

3 time treatments as 1 procedure.

The first treatment is to destroy the fat cells.

The second treatment is to reduce the fat.

The third treatment is to keep the treatment result.

The new technology is used to reduce the fat,not like some machines to reduce the body water, real treat result for reducing the fat,loss weight, keep slim.

**What is Cold laser?**

Cold laser is a pain-free, non-invasive device used for reducing fat by shining the emitted light (laser) over the skin. Cold laser is a revolutionary new fat reduction technology.

A smart alternative to surgery, Cold laser uses a low-level laser to target the stubborn spots that don't respond to regular diet and exercise.

Cold laser gives fast results, with improvements showing immediately after the first treatment and over the course of treatments. Cold laser is ideal for someone who has made a decision to improve their overall diet and exercise regime, or has struggled to get results from a particular area despite regular exercise and healthy eating.

Cold laser is a spot fat reduction technique that help lose weight in the areas that are traditionally much harder to target.

**How it Works?**

Cold laser targets fat in specific problem areas. Laser pads are positioned on the target area such as upper arms, abdomen or thighs. This is a big advantage over diet and exercise which can reduce overall body fat but not shape individual areas.

**Treatment**

A standard treatment requires laser diodes to be placed over the treated area, with an additional two diodes placed over the lymphatic region.

With the pads in place the client is then free to enjoy any additional services or treatments on offer during their 20 minute treatment.

**Does it Hurt?**

Cold laser fat reduction is a low level laser also sometimes known as a cold laser because most people feel very little when it is used on them. Cold laser is a comfortable relaxing treatment where the typical sensation felt is just a slight warmth where the pads are in contact with your skin.

## SH652 Cavitation RF Diode Laser Liposuction Machine



### Specifications:

Wavelength: 650nm  
 Lader Handles: 4 big pads, 2 smaller pads,  
 one Cavitation handle,  
 Cavitation Frequency: 40Khz  
 One six-polar RF Handle  
 RF Frequency :0.5-10KHz  
 RF Energy: 50J  
 Single output: 50nw  
 Number of Lights: 38  
 Screen: Color touch screen

Packing Size: 66\*50cm\*35cm

N.W: 15KG

G.W: 18KG

**Work Principle:**

The Laser-lipolysis emits low levels of laser energy, which creates a chemical signal in the fat cells, breaking down the stored triglycerides into free fatty acids and glycerol and releasing them through channels in the cell membranes. The fatty acids and glycerol are then transported around the body to the tissues that will use them during metabolism to create energy. This process of fatty acid release is a natural response of the body when the body needs to use stored energy reserves, thus Laser-lipolysis is not creating any unnatural reaction in the body nor does it affect or damage any surrounding structures such as skin, blood vessels and peripheral nerves. A period of exercise post treatment will ensure the complete metabolism and thus elimination from the body of the freed fatty acids.

**Accessories list**

1. Host 1set
2. Power supply line 1pc
3. User manual 1copy
4. Laser 2X work head 2pc
5. Laser 4X work head 4pc
6. 40K cavitation head 1pc
7. RF head 1pc
8. Bind ties 1set
9. Accessory frame 1pc
10. Hexagonal screws 4pcs
11. Hexagon wrench 1pc

## SH653 Newest 650nm Diode Laser RF Cavitation Equipment



### Specifications:

Wavelength: 650nm

Handles: 4 big pads, 2 smaller pads,

one Cavitation handle,

One six-polar RF Handle

Single output: 50nw

Number of Lights: 38

Screen: Color touch screen

Size: 52cm x 48cm x 32cm

N.W: 8.5KG

G.W: 11.5KG

**Work Principle:**

The Laser-lipolysis emits low levels of laser energy, which creates a chemical signal in the fat cells, breaking down the stored triglycerides into free fatty acids and glycerol and releasing them through channels in the cell membranes. The fatty acids and glycerol are then transported around the body to the tissues that will use them during metabolism to create energy. This process of fatty acid release is a natural response of the body when the body needs to use stored energy reserves, thus Laser-lipolysis is not creating any unnatural reaction in the body nor does it affect or damage any surrounding structures such as skin, blood vessels and peripheral nerves. A period of exercise post treatment will ensure the complete metabolism and thus elimination from the body of the freed fatty acids.

**Accessories list**

- 1.Host 1set
- 2.Power supply line 1pc
- 3.User manual 1copy
- 4.Laser 2X work head 2pc
- 5.Laser 4X work head 4pc
- 6.40K cavitation head 1pc
- 7.RF head 1pc
- 8.Bind ties 1set
- 9.Accessory frame 1pc
- 10.Hexagonal screws 4pcs
- 11.Hexagon wrench 1pc

## SH654 Latest Cavitation RF Diode Laser Lipolaser



### Specifications:

Wavelength: 650nm

Handles: 4 big pads, 2 smaller pads,

one Cavitation handle,

One six-polar RF Handle

Single output: 50nw

Number of Lights: 38

Screen: Color touch screen

Size: 52cm x 48cm x 32cm

N.W: 8.5KG

G.W: 11.5KG

**Work Principle:**

The Laser-lipolysis emits low levels of laser energy, which creates a chemical signal in the fat cells, breaking down the stored triglycerides into free fatty acids and glycerol and releasing them through channels in the cell membranes. The fatty acids and glycerol are then transported around the body to the tissues that will use them during metabolism to create energy. This process of fatty acid release is a natural response of the body when the body needs to use stored energy reserves, thus Laser-lipolysis is not creating any unnatural reaction in the body nor does it affect or damage any surrounding structures such as skin, blood vessels and peripheral nerves. A period of exercise post treatment will ensure the complete metabolism and thus elimination from the body of the freed fatty acids.

**Accessories list**

- 1.Host 1set
- 2.Power supply line 1pc
- 3.User manual 1copy
- 4.Laser 2X work head 2pc
- 5.Laser 4X work head 4pc
- 6.40K cavitation head 1pc
- 7.RF head 1pc
- 8.Bind ties 1set
- 9.Accessory frame 1pc
- 10.Hexagonal screws 4pcs
- 11.Hexagon wrench 1pc



## SH911 Vacuum Roller Bipolar RF 940nm Infrared Laser



### Treatment Principle

1. Infrared laser lowers skin impedance by heating skin and RF energy penetrates deeply into connective tissue.
2. The synergistic combination of infrared laser and conducted RF energy increases the oxygen intracellular diffusion by heating skin.
3. The technology that vacuum folds skin makes RF energy penetrate a specific folded skin, significantly improving effect and safety, even for the upper eyelid area treatment.
4. Vacuum plus specially-designed rollers manipulate leads even RF penetration (5-15mm). Meanwhile, vacuum and roller mechanical tissue manipulation nips and stretch fibrillar connective tissue, effectively break down subcutaneous fat as well as extruded capillary

vessel, increase lymphatic drainage, promote metabolism and reduce or shrink the size of the actual fat chamber, thereby greatly improving body shaping effect.

**Specification:**

- 1) screen,
  1. instrument screen, 10.4"FT
  2. The screen for hand piece,
    - The screen for hand piece 1,2.4"
    - The screen for hand piece 2,1.9"
- 2) Working mode, pulse sequence  
Pulse spread, 0.5s-7.5s
- 3) Negative pressure,
  1. absolute value, 80kPa -10kPa(60.8Hg -7.6Hg);
  2. the relative value, 20kPa -90kPa(15.2Hg -68.4Hg);
- 4) Roller speed, 0-36rpm
- 5) Roller working Pattern, 2 kinds
- 6) safety checking, On line
- 7) RF frequency, 1MHz
- 8) RF energy density, Max60J/cm
- 9) Laser wavelength, 940nm
- 10) Laser power, Max 20W
- 11) Hand piece quantity, 4
- 12) Therapy square measure, 4mm7mm,8mmx25mm,30mmx50mm,40mmx60mm
- 13) IC card System
- 14) Input power , 750W
- 15) Power supply, AC220V/AC110V,50Hz/60Hz
- 16) Net weight, 50kg
- 17) Size, 65cm\*56cm\*150cm

**Warranty , 1 year ,**

1. Within 1 year, if there is some accessory going wrong, we will send you the accessory at once, free of charge.
2. Within 1 year, if the machine is of any big fails, we will send you a new machine.

# Vacuum Cavitation RF

## SH822 Portable Body RF Face Cavitation Vacuum Roller



### **SPECIFICATION**

Type: RF+Cavitation+Vacuum Roller

Handle: Cavitation+Vacuum RF+Body RF+Face RF

Display: LCD 7 inches touching screen

Cavitation Power: 150W

Cavitation Frequency: 40Khz

Pressure for vacuum: 0.1 Mpa

Vaccum Intensity : 50L/min

RF Frequency :0.5-10KHz

RF Energy: 50J

Operative mode: continuous and pulsed

Power supply: 220V/110V, 50Hz/60Hz

G.W: 26kg

MEAS:75×58cm×49cm

#### **WORK THEORY**

The equipment takes the strong ultrasonic heads to release 40,000Hz powerful ultrasonic, speedy vibrates fatness cells, produces countless vacuum air pocket in and out of cells, mightily impacts the fatty cells, makes the cell membrane produce the inner cracking, so as to dissolve the triglyceride to be the glycerin and free fatty acid, gets the body slimming and beauty effects.

#### **Handle 1. Vacuum Roller**

**1.Stereo vacuum negative pressure, mechanical scroll and skin folding & RF frequency**

**2. 40K Cavitation**

**3. Tripolar RF (Body)**

**4.Tripolar RF( Face)**

**5. Obvious effect**

**6. System Selection Interface**

#### **FUNCTION AND THERAPEUTIC THEORY**

1. Mightily cracking the cellulite, remove the fatness.

It uses the strong ultrasonic head to shoot 40,000Hz ultrasonic, after entering the body, will make the fatty cells produce mighty cracking effect and fatty cells friction, effectively exhaust the heat energy, exhaust cells water, let fatty cells shrinking. Moreover, when the ultrasonic vibrates, there will bring mighty impaction among cells, fatty cells will be less, so get effect of fatness removal.

2. Dissolve the fatness, do lymph drainage, skin tightening, enhance skin elasticity.

3. Improve the orange tissues, slim and sculpture the body contour.

Through the electrode energy of the equipment, it produces biology waves to stimulate the body relevant points. With different frequencies and pulses, under many physics electron reciprocity, stimulate the fatness, let the body exercise, so exhaust the heat energy and fatness again, get the effect of sculpturing the body.

**Warranty : 1 year**

**Shipping way, by DHL, UPS, EMS, Aircargo**

## **GS6.8 Body Analyse with Vacuum Cavitation Slimming Instrument**



**Brief introduction,**

Fast skin-tighten instrument is the perfect system of body shaping, which is auto-fat composition scientific analyse the shape body of vacuum, ultrasound fat lipolysis, bipolar bio-Firming and skin tighten etc.

**Specification:**

Power: <250W,

Cold light: RED GREEN BLUE AND PURPLE

Test age: 10-90 years old;

Printer: USB,

Parallel micro-print

Ultrasound Frequency: 1MHZ

Packing. 68CM\*53CM\*146.5CM,

GW:40KG, NW:31.1KG;

Voltage: AC110, AC220

**Function:**

1. Use the body's eight biological resistance, 5-50KHZ scientific way of calculation, accurate analyse the fat of human body fat, work out the suitable shaping body program for the guests.
2. The meticulous vacuum body slimming system , can shift fat and eliminate fat, stimulate blood circulation, restructuring the organization of fiber;
3. One million times of the ultra-highspeed vibration, melting the body of surplus fat, break down the fat of accumulation;
4. Micro-biological current promote lymph flow, tighten relax skin, improve the edema, accelerate skin regenerate.
5. Touch screen design, patent free-class regulator, the super pressure,
6. Intelligent IC control system, the original pump of America;
7. Vacuum working head for different parts, Fat-testing rod, 1MHz supersonic fat-dissolving head and also the BIO polar skin tighten rod.

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

## GS8.1 Vacuum Cavitation Slimming beauty Equipment



**Product name:**

Vacuum Cavitation Slimming beauty Equipment

**Technical specification**

- 1.Handle: 7 handles
- 2.Output Power :0.5-25W/cm2
- 3.Frequency:100KHZ/80KHZ/40KHZ
- 4.Vacuum pressure:0.1 Mpa
- 5.Power:0.5-25w/cm2
- 6.Gross weight:70kg
- 7.packing, 40\*57\*138cm

**The scope of treatment**

- 1.promote tissue metabolism.
- 2.repel the cellulite.
- 3.tighten the skin
- 4.strengthen the skin elasticity
- 5.fat dissolving
- 6.body shaping
- 7.Net weight:40KG
- 8.Dimension:126x55x58cm
- 9.Voltage: AC110V,AC220V

**Theory of operating**

vacuum supersonic fat-explosion machine can improve tissue metabolism, repel the cellulite, tighten the skin, strengthen the skin elasticity and so on. its effect will last for a long time. At the meanwhile, it can enhance the muscle elasticity and improve the cell metabolism.

supersonic wave creates an effect like internal explosion. it means that the expansion and compression of wave produce a lot of tiny space inside the liquid. the liquid and vapour exist in these tiny space. in the compression cycle, supersonic wave impose a positive pressure on liquid molecule, whereas a negative pressure will be created through the expansion cycle. there is a cohesion effect inside the liquid or biological tissue. inside the low-density lipocyte, the cohesive effects are comparatively weak. A negative pressure, which is relative low, can create tissue gaps easily. it is called as cavitations in physics. the explosion effects inside and outside the cells will increase the molecular movements, so as to achieve the high energy state. then lead to the lipocyte cracking

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.



## GS8.2 E Portable Cavitation Machine



### The description of production

Cavitation Body Slimming machine can promote tissue metabolism, repel the cellulite, tighten the skin, strengthen the skin elasticity and so on. Its effect will last for a long time. Meanwhile, it can enhance the elasticity of muscle and improve the metabolism of cells.

Supersonic wave create an effect like internal explosion. It means that the expansion and compression of the wave produce a lot of tiny space inside the liquid. The gas and vapour exist in these tiny space. In the compression cycle, supersonic wave imposes a positive pressure on liquid molecule, whereas a negative pressure will be created through the expansion cycle. There is a cohesion effect inside the liquid and biological tissues. Inside the low-density lipocyte, the cohesive effects are comparatively weak. A negative pressure, which is relative low can create tissue gaps easily. It is called as "cavitations" in physics. The

explosion effects inside and outside the cells will increase the molecular movements, so as to achieve an high energy state, and then lead to the lipocyte cracking.

**Specifications:**

- 1.Handle :3 handles
- 2.Out put power: 5-25W/cm<sup>2</sup>
- 3.Frequency: 40KHZ/1MHZ/1MHZ
- 4.Controller: liquid crystal button screen
- 5.Gross weight: 15kgs
- 6.Voltage: AC110V,AC220V

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 years warranty, Within warranty, we offer free repairing service and provide replacement parts for free.



## M9 Super Body Sculptor



### Main Theory

Love of the crystal is an advanced instrument of cosmetology and bodyfitting in the 21st century, because it fuses the medical cosmetology electricity face-lifting, the RFpen of making the skin supple, the head of medical electricity face-lifting, the ice wave head, exploiting the fat head, the radio frequency head of losing weight, whose function is the following: removing the dark pouches, the pouch and the wrinkle of eyes; eliminating the wrinkle of face, promoting, whitening and making the skin tender, wiping off the scar; losing weight with high efficiency, removing the striae gravidarum and tightly tightening the flesh. Because of its merits such as its effect is quick and lasting, non-recoil and operation is simple, it is popular with beautician and people who love beauty.

### Specification

1. Input voltage: AC 110V/220V 50/60HZ
2. Input power: Less than 650W

- 3.Output frequency:0.2-10Khz
- 4.Output RF:5MHZ
- 5.Output pressure:0-250Kpa
- 6.Size:96\*121\*154CM

**Features**

- 1 video detection system(proprietary technology)
- 2 video monitoring system
- 3 60°Fat rotating system(proprietary technology)
- 4 5M multipolar RF thermal slimming system
- 5 Vacuum & photon move fat system
- 6 35K ultrasonic cavitation system
- 7 heart rate test technology
- 8 Color touch screen with multimedia, display and control the procedure
- 9 USB data transmit

The design and the current is fully comply with European electrical safety standards, absolutely safe to human body, the built-in anti-electromagnetic interference devices, can effectively protect against electrical and magnetic interference, maintain the current stability, to ensure that procedures effective and stable

**Five system:**

- Vacuum optional fiber fat-sucking system
- RF 450K high-frequency fat-dissolving system
- Energy and fat rotary system
- Processing video monitoring system
- Body and face shaping,weight loss

**Function**

slimming equipment lose weight ,vacuum strong can promote tissue metabolism,repel the cellulue.tighrten the skin ,strengthen the skin elastic and so on ,its effect will last for a long time.at the same time ,it can enhance the muscle elasticity and improve the cell metabolism with rolers in it ,can move different directions and different speed

- 1.Replacing and sumouning the surgery of attracting the fat;
- 2.Tightening the flaccid skin;
- 3.Promoting the speed of organism metabolism and accelerating the waste and the moisture extraction of the body;
- 4.Reparing striae gravidarum;
- 5.Relaxing themuscies,relieving the muscle spasm and alleviating the ache symptomof muscles;
- 6.Losing weight of the whole body.

## GS9.2 Vacuum Cavitation Slimming Beauty Machine



### Technical specification

1. Display: LCD 10.4 inches touching screen
2. Frequency: 40Khz/80Khz/100Khz
3. Output power: 0.5-50W/cm<sup>2</sup>
4. Handle: 7 handles: 1 M80 + 1 M50 + 4X + 1 vacuum
5. Output energy: M80 is 20-50 W/cm<sup>2</sup>  
M50 and 4X is 5-15W/cm<sup>2</sup>  
Vacuum is 0.5-5W/cm<sup>2</sup>
6. Operative mode: continuous and pulsed
7. Pressure for vacuum: 0.1mkPa.
8. Shell material: stainless steel
9. Power supply: 220V/110V 50Hz/60Hz
10. Fuse: 2x5.15 AT

11. Temperature: 10-45°C
12. Humidity: 0-85%
13. Weight: 60Kg
14. Packing, 38\*68\*128cm

**General introduction**

Ultrasound easy to produce a liquid implosion effect, Namely, wave expansion and compression form a large number of micro-gap in the liquid, which is full of gas and steam, the strong sound waves have positive effects on the liquid molecules at compression cycle, while it has vacuum effect at expansion cycle. There is cohesion among liquid and biological tissues, molecular bonding is weak in low density fat cells, and the low vacuum caused by strong sound waves can generate organized gaps, in physics known as the "cavitations". And the implosion caused by micro gaps inside and outside cells will promote molecular motion, makes a high energy level, which finally led to break of fat cells. Vacuum & ultrasound slimming machine applies the cavitation theory, uses the 40kHz~1MHz ultrasound to make intense vibration, which causes an inner explosion in the fat tissues. thus get rid of lipocyte. It can effectively break down stubborn and deep fat, decrease lipocyte quantity and its size. The diathermal and vacuum suction can promote lipocyte and tissues liquid flow, which help drainage of fatty acid, waste and toxin through lymph system. The vacuum function will promote liquid flow of lipocyte and fat tissues, then speed up excretion of toxins and waste, achieve effect of fat elimination, slimming, skin tightening and enhancing muscle flexibility.

**Application Range**

1. promote tissue metabolism.
2. repel the cellulite.
3. tighten the skin
4. strengthen the skin elasticity
5. fat dissolving
6. body shaping

## GS9.3 Vacuum Cavitation Goethe-shi Electricity Current Beauty Equipment



### Specifications:

- 1.Goethe-shi electricity current
- 2.40KHz ultrasound, 50KHz ultrasound
- 3.Cavitation heads+Vacuum head
- 4.Output Power: 0.5-25W/cm<sup>2</sup>
- 5.Handles: Four heads+Electrode Plate , one 40K cavitation head, one 50K cavitation head
- 6.2 different sizes of vacuum heads for small and big areas
- 7.Screen: Color touch screen
- 8.Host Power: <150W

- 9.Vacuum: 0.1Mpa
- 10.Frequency 1 : 1-250Hz
- 11.Frequency 2 : 40KHz
- 12.Frequency 3 : 50KHz
- 13.Dimensions: 36\*53\*96cm
- 14.Net Wight: 35KG

**Applications:**

Vacuum and cavitation system,

Vacuum strong sound wave fat system,can promote tissue metabolism,repel the cellulite,tighten the skin,elastic and so on.Its effect will last for long time.


**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.





## RU+2 Tripolar RF Cavitation Equipment

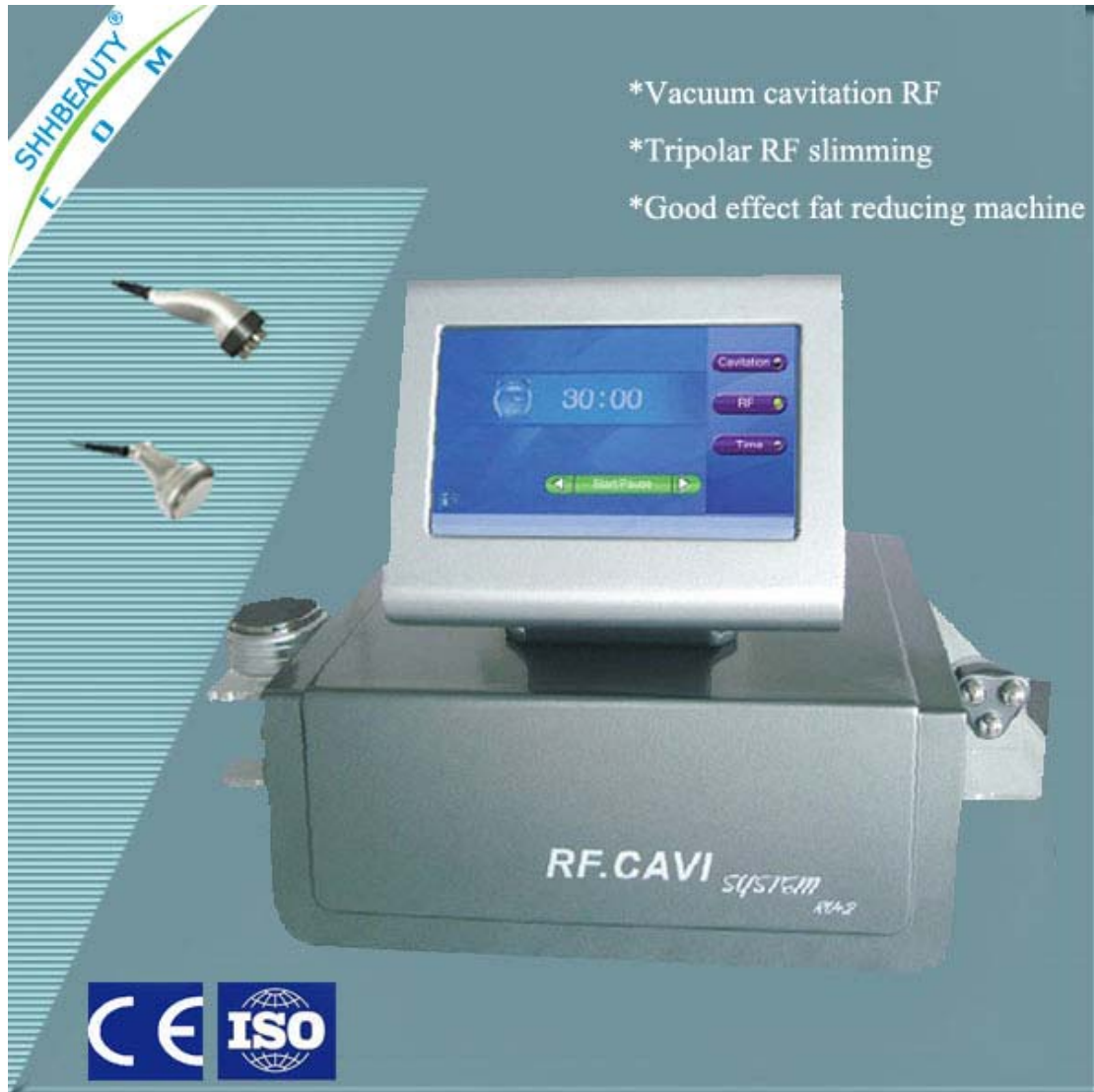

### SPECIFICATION

Output Power : 110W

big 10.4 color touch screen

Frequency:40KHZ 5W-50W deepest fat

tripolar RF 5MHz Various poles

Voltage :110-230V/50/60hz

### TECHNICAL ADVANTAGES

Adopt the top ultrasonic liposuction technique in the world.

Suitable for all kinds of skin.

Comfortable, painless, woundless during the treatment.

Easy operation, easy study.

Non- exhaustion, low cost investment and quick high returns..

### **FEATURES**

Free operation, . non-anaesthetic.

No ruggedness.

No bleeding, tumidness and bruise.

No side effects, good effects, no rebounding phenomenon.

Woundless, will not influence the normal working and living.

### **WORK THEORY**

The equipment takes the strong ultrasonic heads to release 40,000Hz powerful ultrasonic, speedy vibrates fatness cells, produces countless vacuum air pocket in and out of cells, mightily impacts the fatty cells, makes the cell membrane produce the inner cracking, so as to dissolve the triglyceride to be the glycerin and free fatty acid, gets the body slimming and beauty effects.

### **FUNCTION AND THERAPEUTIC THEORY**

1. Mightily cracking the cellulite, remove the fatness.

It uses the strong ultrasonic head to shoot 40,000Hz ultrasonic, after entering the body, will make the fatty cells produce mighty cracking effect and fatty cells friction, effectively exhaust the heat energy, exhaust cells water, let fatty cells shrinking. Moreover, when the ultrasonic vibrates, there will bring mighty impaction among cells, fatty cells will be less, so get effect of fatness removal.

2. Dissolve the fatness, do lymph drainage, skin tightening, enhance skin elasticity.

3. Improve the orange tissues, slim and sculpture the body contour. Through the electrode energy of the equipment, it produces biology waves to stimulate the body relevant points. With different frequencies and pulses, under many physics electron reciprocity, stimulate the fatness, let the body exercise, so exhaust the heat energy and fatness again, get the effect of sculpturing the body.

### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.

2. The professional and innovative engineers are the pioneers of optoelectronics in China.

3. ISO9001 proved Quality control and test system, each machine is in the best quality control.

4. CE proved machines, the key for you to legally use and sell the machine.

5. Spare parts self produced, production cost is in control.

6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

## RU+5 Multipolar RF Vacuum Cavitation Slimming Machine



### Specification

Output Strong sound waves: 40KHz/1MHz;

Output RF: 5MHz

Input Voltage: AC110V/220V 50-60HZ

Power: <500w

NW.: 14KGS;

GW.: 20KGS

Meas: 53CM\*70CM\*65cm

**RU+5 features:**

slimming equipment lose weight ,vacuum strong can promote tissue metabolism,repel the cellulite.tighten the skin ,strengthen the skin elastic and so on ,its effect will last for a long time.at the same time ,it can enhance the muscle elasticity and improve the cell metabolism with rollers in it ,can move different directions and different speed

- 1.Replacing and sumouning the surgery of attracting the fat;
- 2.Tightening the flaccid skin;
- 3.Promoting the speed of organism metabolism and accelerating the waste and the moisture extraction of the body;
- 4.Reparing striae gravidarum;
- 5.Relaxing themuscies,relieving the muscle spasm and alleviating the ache symptomof muscles;
- 6.Losing weight of the whole body.

**RF vacuum cavitation slimming machine main Theory:**

Love of the crystal is an advanced instrument of cosmetology and bodyfitting in the 21st century,because it fuses the medical cosmetology electricity face-lifting,the RFpenof making the skin supple,the head of medical electricity face-lifting,the ice wave head,exploing the fat head,the radio frequency head of losing weight,whose function is the following:removing the dark pouches,the pouch and the wrinkle of eyes;eliminating the wrinkle of face,promoting,whitening and making the skin tender,wiping off the scar;losing weight with high efficiency,removing the striae gravidarum and tightly tightening the flesh.Because of its merits such as its effect is quick and lasting,non-recoil and operation is simple,it is popular with beautician and people who love beauty.

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

# RU+6 Cavitation Multipolar RF Tripolar RF Bipolar RF Slimming Machine



## Accessories



Bipolar RF work head



Tripolar RF work head



Hexapolar RF laser work head



40K cavitation explode work head



Service brake

**Technical parameter**

AC220v/110v 50-60Hz  
power: <500w  
output ultrasonic:40KHz/1MHz  
output RF:5MHz  
vacuum pressure:100KPa

**Functions**

- 3 polar RF handle for face treatment.
- 6 polar handle for body treatment.
- 2 polar working head for eyes
- M70 cavitation handle for fat burning

**Working theory**

it creates an effect like internal Explosion. It means that the expansion and compression of wave produce a lot of tiny space inside the liquid. The gas and vapour exist in these tiny space. In the compression cycle, supersonic wave impose a positive pressure on liquid molecule, whereas a negative pressure will be created through the expansion cycle. There is a cohesion effect inside the liquid or biological tissues. Inside the low-density lipocyte, the cohesive effects are comparatively weak. A negative pressure, which is relative low, can create tissue gaps easily. It is called as Cavitations in physics. The explosion effects inside and outside the cells will increase the molecular movements, so as to achieve an high energy state, then lead to the lipocyte cracking.

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

## SH2.8 Mini Cavitation Vacuum Bipolar RF Multi RF RF Slimming Machine



### SPECIFICATION

Type: RF+Cavitation+Vacuum

Handle: Cavitation+Vacuum RF+Body RF+Face RF

Display: LCD 7 inches touching screen

Cavitation Power: 100W

Cavitation Frequency: 40Khz

Pressure for vacuum: 0.1 Mpa

RF Frequency: 1MHZ

RF Energy: 50J

LED Light: red-630nm, blue-415nm

Operative mode: continuous and pulsed

Size: 54\*48\*39cm

Power supply: 220V/110V 50Hz/60Hz

G.W: 22kg

#### **TECHNICAL ADVANTAGES**

Adopt the top ultrasonic liposuction technique in the world.

Suitable for all kinds of skin.

Comfortable, painless, woundless during the treatment.

Easy operation, easy study.

Non- exhaustion, low cost investment and quick high returns..

#### **WORK THEORY**

The equipment takes the strong ultrasonic heads to release 40,000Hz powerful ultrasonic, speedy vibrates fatness cells, produces countless vacuum air pocket in and out of cells, mightily impacts the fatty cells, makes the cell membrane produce the inner cracking, so as to dissolve the triglyceride to be the glycerin and free fatty acid, gets the body slimming and beauty effects.

#### **FUNCTION AND THERAPEUTIC THEORY**

1. Mightily cracking the cellulite, remove the fatness.

It uses the strong ultrasonic head to shoot 40,000Hz ultrasonic, after entering the body, will make the fatty cells produce mighty cracking effect and fatty cells friction, effectively exhaust the heat energy, exhaust cells water, let fatty cells shrinking. Moreover, when the ultrasonic vibrates, there will bring mighty impaction among cells, fatty cells will be less, so get effect of fatness removal.

2. Dissolve the fatness, do lymph drainage, skin tightening, enhance skin elasticity.

3. Improve the orange tissues, slim and sculpture the body contour.

Through the electrode energy of the equipment, it produces biology waves to stimulate the body relevant points. With different frequencies and pulses, under many physics electron reciprocity, stimulate the fatness, let the body exercise, so exhaust the heat energy and fatness again, get the effect of sculpturing the body.



**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

**SH2.9 Cavitation+Vacuum RF+Body RF+Face  
RF+BIO**



## **SPECIFICATION**

Type: RF+Cavitation+Vacuum+BIO

Handle: Cavitation+Vacuum RF+Body RF+Face RF+BIO

Display: LCD 7 inches touching screen

Cavitation Power: 100W

Cavitation Frequency: 40Khz

Pressure for vacuum: 0.1 Mpa

RF Frequency: 1MHZ

RF Energy: 50J


LED Light: red-630nm, bule-415nm

Operative mode: continuous and pulsed

Size: 53\*48\*49cm

Power supply: 220V/110V 50Hz/60Hz

G.W: 23kg

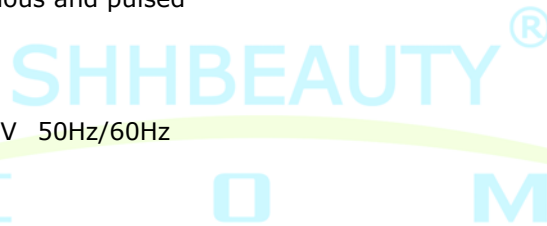

## **TECHNICAL ADVANTAGES**

Adopt the top ultrasonic liposuction technique in the world.

Suitable for all kinds of skin.

Comfortable, painless, woundless during the treatemnt.

Easy operation, easy study.

Non- exhaustion, low cost investment and quick high returns..

## **WORK THEORY**

The equipment takes the strong ultrasonic heads to release 40,000Hz powerful ultrasonic, speedy vibrates fatiness cells, produces countless vacuum air pocket in and out of cells, mightily impacts the fatty cells, makes the cell membrane produce the inner cracking, so as to dissolve the triglyceride to be the glycirin and free fatty acid, gets the body slimming and beauty effects.

## **FUNCTION AND THERAPEUTIC THEORY**

1. Mightily cracking the cellulite, remove the fatness.

It uses the strong ultrasonic head to shoot 40,000Hz ultrasonic, after entering the body, will make the fatty cells produce mighty cracking effect and fatty cells friction, effectively exhaust the heat energy, exhaust cells water, let fatty cells shrinking. Moreover, when the ultrasonic vibrates, there will bring mighty impaction among cells, fatty cells will be less, so get effect of fatness removal.

2. Dissolve the fatness, do lymph drainage, skin tightening, enhance skin elasticity.

3. Improve the orange tissues, slim and sculpture the body contour.

Through the electrode energy of the equipment, it produces biology waves to stimulate the body relevant points. With different frequencies and pulses, under many physics electron reciprocity, stimulate the fatness, let the body exercise, so exhaust the heat energy and fatness again, get the effect of sculpturing the body.

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.

2. The professional and innovative engineers are the pioneers of optoelectronics in China.

3. ISO9001 proved Quality control and test system, each machine is in the best quality control.

4. CE proved machines, the key for you to legally use and sell the machine.

5. Spare parts self produced, production cost is in control.

6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.

## SH2.11 Cavitation Tripolar RF Multipolar RF Machine



### Technical parameters:

- 1.RF frequency: 1-10MHz
- 2.RF energy density: Max 60J/cm<sup>3</sup>
3. focused blast frequency FAT: 35KHz-45KHz
4. power: 50W
5. power supply: AC230V ± 10%,50Hz ± 1Hz/AC110V ± 10%,60Hz ± 1Hz
6. weight: 8kg
- 7.CTN dimension: 42\*38\*24CM

### The device profile

Using "void" principle, 40kHz strong sound wave frequency emission from strong vibration, blast effect in the adipose tissue, fat natural blast because the subject does not mean inside and outside of the cell membrane, effectively

break down stubborn fat and explode deep fat, reducing fat cell size and number of fat cells. Strong sound wave heating effect of strengthening the flow of liquid between the fat cells and tissues, will blast fat acids, waste and toxins from the lymph system is eliminated from the body, to reduce fat, thin body, skin tightening, enhance muscle elasticity effects.

## Elight RF Nd Yag Laser

### SH1.0 E-light RF Nd Yag Laser Multi-Function Beauty Equipment



### Specifications

- 1.Wavelength:430nm/480nm/530nm/580nm/590nm/640nm/690nm(optonal)
- 2.Operating interface:8.4 inch color touch screen Laser:5.1 inch Display screen
- 3.Frequency:1.2.3HZ
- 4.Elght energy:0-50J/CM2(IPL),0-50J/CM2(RF)
- 5.Laser energy level:1-100 level
- 6.Output power:Elght:800W Laser:350W
- 7.Laser treatment head:1064nm,532nm,1320nm,moppet head
- 8.Cooling system:Semiconductor+wind cooling+water cooling
- 9.Input power:AC 220V(110V)/15A 50HZ(60HZ)
- 10.Machine size:40\*53\*123cm
- 11.Package size:51\*70\*134cm

### Application

- 1.Permanent hair removal by E-light technology
- 2,All kinds of spots removal, freckle removal, pigment removal
- 3.Skin lifting, skin tightening, skin rejuvenation, skin whitening, wrinkle removal, reshape facial outlines, contract large pores, rebuild body figure
- 4.Acne treatment, Angiotelectas treatment
- 5.All colors of Tattoo removal, eyebrow line, lip line pigments removal
- 6.Weight loss, breast enlarging.

### Features

1. A good combination of five advanced technologies: IPL+E-light+ Bi-polar RF pigment removal+ Laser tattoo removal +Soft Skin Black Moppet, it is the KING of all the beauty machines, which could solve nearly all the skin problems.
2. Double operation screen: 8.4 inch TFT True colorful touch screen+5.1" LCD screen, which could work together, easy to operate.
3. Fashionable style, very popular to customers, easy to get the invest back.
4. CPC connector of instant plug and play, very easy to install and operate.
5. Modular Design of inner structure, which makes it work more stable.

### Company advantages:

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

## SH1.5 E-light IPL RF Nd Yag Laser



### Specifications,

- 1.Wavelength:430nm/480nm/530nm/580nm/590nm/640nm/690nm(optonal)
- 2.Operating interface:8.4 inch color touch screen Laser:5.1 inch Display screen
- 3.Frequency:1.2.3HZ
- 4.Elight energy:0-50J/CM2(IPL),0-50J/CM2(RF)
- 5.Laser energy level:1-100 level
- 6.Output power:Elight:800W Laser:350W
- 7.Laser treatment head:1064nm,532nm,1320nm,moppet head
- 8.Cooling system:Semiconductor+wind cooling+water cooling
- 9.Input power:AC 220V(110V)/15A 50HZ(60HZ)
- 10.Machine size:40\*53\*123cm
- 11.Package size:51\*70\*134cm

### **Application**

1. Permanent hair removal by E-light technology
2. All kinds of spots removal, freckle removal, pigment removal
3. Skin lifting, skin tightening, skin rejuvenation, skin whitening, wrinkle removal, reshape facial outlines, contract large pores, rebuild body figure
4. Acne treatment, Angiotelectas treatment
5. All colors of Tattoo removal, eyebrow line, lip line pigments removal
6. Weight loss, breast enlarging.

### **Features**

1. A good combination of five advanced technologies: IPL+E-light+ Bi-polar RF pigment removal+ Laser tattoo removal +Soft Skin Black Moppet, it is the KING of all the beauty machines, which could solve nearly all the skin problems.
2. screen: 10.4 inch TFT True colorful touch screen
3. Fashionable style, very popular to customers, easy to get the invest back.
4. CPC connector of instant plug and play, very easy to install and operate.
5. Modular Design of inner structure, which makes it work more stable.

### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.



# SH1.6 E-light IPL RF Nd Yag Laser Multi-functional Beauty Machine



## Specifications,

1.Wavelength:430nm/480nm/530nm/580nm/590nm/640nm/690nm; 2.Operating interface:10.4h color touch screen; 3.Frequency:1.2.3HZ;; 4.Elight energy:0-50J/CM2(IPL),0-50J/CM2(RF); 5.Laser energy level:1-100 level; 6.Output power:Elight:800W Laser:350W; 7.Laser treatment head:1064nm,532nm,1320nm,moppet head; 8.Cooling system:Semiconductor+wind cooling+water cooling; 9.Input power:AC 220V(110V)/15A 50HZ(60HZ); 10.Machine size:40\*53\*123cm; 11.Package size:60\*68\*138cm

## Application

- 1.Permanent hair removal by E-light technology
- 2,All kinds of spots removal, freckle removal, pigment removal

- 3.Skin lifting, skin tightening, skin rejuvenation, skin whitening, wrinkle removal, reshape facial outlines, contract large pores, rebuild body figure
- 4.Acne treatment, Angiotelectas treatment
- 5.All colors of Tattoo removal, eyebrow line, lip line pigments removal
- 6.Weight loss, breast enlarging.

**Features**

1. A good combination of five advanced technologies: IPL+E-lightI+ Bi-polar RF pigment removal+ Laser tattoo removal +Soft Skin Black Moppet, it is the KING of all the beauty machines, which could solve nearly all the skin problems.
2. screen: 10.4 inch TFT True colorful touch screen
2. Fashionable style, very popular to customers, easy to get the invest back.
4. CPC connector of instant plug and play, very easy to install and operate.
5. Modular Design of inner structure, which makes it work more stable.

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

## SH880 DIY Beauty Equipment



### Specification:

Voltage: 220V/110V

Frequency 50/60 Hz

Output: 150W

Gross Weight: 60kg

Net Weight: 35kg

Size: 145\*74\*74cm

### Working Theory:

The equipment takes the strong ultrasonic heads to release 40,000Hz or 1 MHz powerful ultrasonic, speedy vibrates fattiness cells, produces countless vacuum air pocket in and out of cells, mightily impacts the fatty cells, makes the cell membrane produce the inner cracking, so as to dissolve the triglyceride to be the glycerin and free fatty acid, then under the function of cavitations, to exhaust them out through the liver sausage circulation, at last, uses the

Ultrasonic Locating plate to do skin tightening, gets the body slimming and beauty effects

**Function**

- 1.Tighten the flabby skin.
- 2.Increase the bodys metabolic rate; help the body excrete waste and moisture quickly.
- 3.Repair striate gravid arum.
- 4.Relax the muscles, lift the muscle cramps, and ease muscle pain and so on.
- 5.Lose weight.
- 6.Effectively cure the orange peel skin on buttock or thigh, solve the problem of the flabby abdominal skin after procreating or liposuction

**Technical Advantages**

- 1.Adopt the top ultrasonic liposuction technique in the world.
- 2.Suitable for all kinds of skin.
- 3.Comfortable, painless, woundless during the treatment.
- 4.Easy operation, easy study.
- 5.Non- exhaustion, low cost investment and quick high returns.

**Optional functions(Option to add):**

- M201:Ultrasonic
- M202:High frequency
- M203:Skin scrubber
- M204:Hot & Cold bar
- M205:Galvanic
- M206:Facial brush



# Nd Yag Laser

## SH3.0 Nd Yag Laser Tattoo Removal Beauty Equipment



### Brief introduction

- Q-Switched Nd:YAG laser uses an intense beam of light to significantly lighten or completely remove your tattoo. When the ink particles of your tattoo absorb the light from the laser, they are broken up into tiny fragments.
- These fragments are then absorbed through your body's natural cleansing

mechanism. This cleansing process generally takes a few weeks. Larger ink particles may require multiple treatments.

#### **Functions**

- Remove kinds of colors tattoos. (blue, brown,black, red, coffee etc tattoos)
- Remove pigment spots, birthmark etc.
- Dispel age spot and beauty spots etc.
- Skin care

#### **Specification**

- Wavelength: 1064nm & 532nm,1320nm
- Laser type: sapphire+ ruby switch Q/KTP/YAG laser instrument
- Energy: 1600 mj
- Indication light: infrared ray indication light
- Width of pulse: 10ns
- Frequency: 1-5 HZ
- Spot diameter: 1-8mm
- Cooling manner: wind+ water
- Input power: AC 220V(110V)/ 15A 50HZ(60HZ)
- Machine Size: 50\*40\*110cm
- Package Size: 56\*53\*116cm

#### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China,
3. ISO9001 proved Quality control and test system, each machine is in the best quality control
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control

#### **Training:**

Once you buy our machines,we supply professional training for you. Detailed user manual and online service are training basis. Moreover, we teach step by step about machine usage and treatment operation.Training certification will be provided after training finish.

#### **After sale service:**

1 year warranty

Within warranty, we offer free repairing service and provide replacement parts for free.

1. we are the only one and the first one to supply 1 year warranty
2. we are very confidence about the quality of our products.
3. we try best to make clients work with the machine without any troubles after buying

## SH3.1 Portable Q-Switched Nd Yag Laser Beauty Machine



### Brief introduction

- Q-Switched Nd:YAG laser uses an intense beam of light to significantly lighten or completely remove your tattoo. When the ink particles of your tattoo absorb the light from the laser, they are broken up into tiny fragments.

•These fragments are then absorbed through your body's natural cleansing mechanism. This cleansing process generally takes a few weeks. Larger ink particles may require multiple treatments.

#### **Functions**

- Remove kinds of colors tattoos. (blue, brown,black, red, coffee etc tattoos)
- Remove pigment spots, birthmark etc.
- Dispel age spot and beauty spots etc.
- Skin care

#### **Specification**

- Wavelength: 1064nm & 532nm,1320nm
- Laser type: sapphire+ ruby switch Q/KTP/YAG laser instrument
- Energy: 1600 mj
- Indication light: infrared ray indication light
- Width of pulse: 10ns
- Frequency: 1-5 HZ
- Spot diameter: 1-8mm
- Cooling manner: wind+ water
- Input power: AC 220V(110V)/ 15A 50HZ(60HZ)
- Machine Size: 24\*28\*42cm
- Package Size: 31\*38\*50cm

#### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China,
3. ISO9001 proved Quality control and test system, each machine is in the best quality control
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control

#### **Training:**

Once you buy our machines,we supply professional training for you. Detailed user manual and online service are training basis. Moreover, we teach step by step about machine usage and treatment operation.Training certification will be provided after training finish.

#### **After sale service:**

1 year warranty

Within warranty, we offer free repairing service and provide replacement parts for free.

1. we are the only one and the first one to supply 1 year warranty
2. we are very confidence about the quality of our products.
3. we try best to make clients work with the machine without any troubles after buying.



# E-light IPL RF

## SH4.0 E-light RF Beauty Equipment



### Technical parameters:

- 1.Light source:Intensive pulse light
- 2.IPL energy:10-50J/cm<sup>2</sup>
- 3.IPL frequency:1s,2s,3s,4s
- 4.Pulse space:5-60ms

- 5.Pulse duration: 2-15ms
- 6.IPL power source:1600W
- 7.RF frequency:6Mhz
- 8.RF power source:600W
- 9.headpiece:Tow
- 10.Light spot area: 12\*30mm<sup>2</sup>/8\*400mm<sup>2</sup>
- 11.Cooling system: semiconductor+ water+ wind+ cooling gel
- 12.input power : AC 220V(110V)/ 15A 50HZ(60HZ)
- 13.Machine Size:50\*53\*133cm

### **Application**

- 1.Permanent hair removal by E-light technology
- 2,All kinds of spots removal, freckle removal, pigment removal
- 3.Skin lifting, skin tightening, skin rejuvenation, skin whitening, wrinkle removal, reshape facial outlines, contract large pores, rebuild body figure
- 4.Acne treatment, Angiotelectas treatment
- 5.Weight loss, breast enlarging.

### **Theory**

E-light=IPL+RF+Skin contact cooling system.

It's the new generation of beauty equipment,it combines with the functions of optical energy and electromagnetic waves.It adopts elective absorption of light,which could get different impedance between target tissue and abnormal skin.It could strengthen absorption of electromagnetic of target tissue in the condition of lower light energy,which could greatly reduce painless during treatment.

### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

## SH4.1 Portable E-light IPL RF Beauty Equipment



### Technical parameters:

- 1.Light source:Intensive pulse light
- 2.IPL energy:10-50J/cm<sup>2</sup>
- 3.IPL frequency:1s,2s,3s,4s
- 4.Pulse space:5-60ms
- 5.Pulse duration: 2-15ms
- 6.IPL power source:1600W
- 7.RF frequency:6Mhz
- 8.RF power source:600W
- 9.headpiece:Tow
- 10.Light spot area: 12\*30mm<sup>2</sup>/8\*400mm<sup>2</sup>

- 11.Cooling system: semiconductor+ water+ wind+ cooling gel
- 12input power : AC 220V(110V)/ 15A 50HZ(60HZ)
- 13.Machine Size:60\*45\*35cm

### **Application**

- 1.Permanent hair removal by E-light technology
- 2,All kinds of spots removal, freckle removal, pigment removal
- 3.Skin lifting, skin tightening, skin rejuvenation, skin whitening, wrinkle removal, reshape facial outlines, contract large pores, rebuild body figure
- 4.Acne treatment, Angiotelectas treatment
- 5.Weight loss, breast enlarging.

### **Theory**

E-light=IPL+RF+Skin contact cooling system.

It's the new generation of beauty equipment,it combines with the functions of optical energy and electromagnetic waves.It adopts elective absorption of light,which could get different impedance between target tissue and abnormal skin.It could strengthen absorption of electromagnetic of target tissue in the condition of lower light energy,which could greatly reduce painless during treatment.



### **Company advantages:**

- 1. Specialized in production of medical and aesthetic equipments for many years.
- 2. The professional and innovative engineers are the pioneers of optoelectronics in China.
- 3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
- 4. CE proved machines, the key for you to legally use and sell the machine.
- 5. Spare parts self produced, production cost is in control.
- 6.We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

# 808 Diode Laser

## SH808 808nm Diode Laser Hair Removal Machine



### Specifications

|                     |  |
|---------------------|--|
| Laser type          | diode laser  |
| Wavelength          | 808nm±10nm   |
| Housing             | stacked array diode laser                                |
| Output power        | 600W(max)  |
| Operation mode      | Single pulse, Double pulse                               |
| Pulse               | Pulse width: 10ms-400ms, Pulse repetition rate: 1Hz-10Hz |
| Transmission system | Crystal Windows  |
| Operation interface | Color LCD touch screen                                   |
| Power supply        | 220Vac, 15A, 50Hz  |

|                             |                                     |
|-----------------------------|-------------------------------------|
| Laser Class                 | 4                                   |
| Safety classification       | I Ib                                |
| Cooling                     | Distilled water+wind+ semiconductor |
| FUSE                        | T 250V 15A,                         |
| Water proof level           | IPX1                                |
| Footswitch Waterproof level | IPX8                                |

(temperature is 5°C~40°C, relative humidity is less then 80%,and atmospherisis is 860hPa~1060hPa)

Machine Packing size: 150\*65\*65cm

Gross Weight: 65kgs

### **Theory and the technical information**

Diode laser is a kind of laser with semiconductor as working material. It consists of working material, cavity resonator and power source.

The diode laser for this unit is GaAlAs diode bar, and the wavelength is 808nm. It features impact structure, high efficiency and long lifetime. Generally the beam shall be shamed as the big beam divergence of the laser from the diode. With the egold'S unique pulse power technology and ICE Colling, the hair can be removal painless.

### **Main Features:**

1. Remove area: Lip hair, beard, chest hair, armpit hair, back hair, arm hair, leg hair, and unwanted hair outside the bikini line.
2. Secure high-speed repair damaged skin.
3. Whitening Firming ultra-efficient, ultra-soft skin can improve
4. Effective treatment for dry skin and rough skin
5. Effective treatment for large pores, fine wrinkles and redness.
6. Effective treatment freckles, reduced skin elasticity and dull skin, etc.

## SH810 808nm Diode Laser Hair Removal Equipment



### Technical parameters

|                       |                                |
|-----------------------|--------------------------------|
| Type of laser         | Continuous semiconductor laser |
| wavelength            | 808nm                          |
| Output                | pulse                          |
| Pulse width           | 0-400ms                        |
| Spot size             | 12mm X 20mm                    |
| Handle cooling system | semiconductor cooling          |
| Frequency             | 1-10HZ                         |
| Energy density        | 0-120J/cm <sup>3</sup>         |
| Laser power           | 500W                           |
| Cooling system        | constant temperature cooling   |
| Input power           | 1800W                          |

## **Brief Introduction**

808 semiconductor laser hair removal instrument can in lower energy density mode will be heated to the hair follicle and °C, through the treatment of hand with low energy sliding a period of time (in the low frequency condition) hair follicle stem cells and growth is lose growth activity, so as to achieve the purpose of permanent hair removal.

## **Working Principle**

Semiconductor freezing point depilatory instrument adopts 808 nm wavelength, the wavelength effective penetration depth can achieve target tissue (hair papilla), appropriate pulse duration security target tissue generate enough thermal damage and surrounding tissues almost not affected; A moderate amount of energy density security in the suitable time to provide strong enough energy output enough to damage the target tissue and normal tissue is almost not affected; Suitable skin protection measures to guarantee the target tissue enough damage, and skin almost not affected, so as to ensure the safety of the treatment.

## **Feature and Advantage**

- 1)Frozen environment , High-speed unhairing reduce pain truly minimally invasive.
- 2)Center wavelength of 808nm through the skin effect hair papilla directly and removal efficiency is higher than other hair removal device.
- 3)Use the 1800W power with the 500W semiconductor laser power supply high frequency and high efficiency light.
- 4)The compressor condenser that handle constant temperature cooling without intermittent work, more suitable for professional medical establishments.
- 5)Using the sliding and function handles ,continuous light sliding operation removal without from the site of action.

## **Application**

Permanent Hair removal of all body part : beard, lips,arm ,leg ,armpit hair,bikini line and so on .Improve pigment ,dark color of the mole ,burgundy birthmarks and other symptoms.



## SH811 808nm Diode Laser



### Technical parameters

|                       |                                |
|-----------------------|--------------------------------|
| Type of laser         | Continuous semiconductor laser |
| wavelength            | 808nm                          |
| Output                | pulse                          |
| Pulse width           | 0-400ms                        |
| Spot size             | 12mm X 20mm                    |
| Handle cooling system | semiconductor cooling          |
| Frequency             | 1-10HZ                         |
| Energy density        | 0-120J/cm <sup>3</sup>         |
| Laser power           | 500W                           |
| Cooling system        | constant temperature cooling   |
| Input power           | 1800W                          |

## **Brief Introduction**

808 semiconductor laser hair removal instrument can in lower energy density mode will be heated to the hair follicle and °C, through the treatment of hand with low energy sliding a period of time (in the low frequency condition) hair follicle stem cells and growth is lose growth activity, so as to achieve the purpose of permanent hair removal.

## **Working Principle**

Semiconductor freezing point depilatory instrument adopts 808 nm wavelength, the wavelength effective penetration depth can achieve target tissue (hair papilla), appropriate pulse duration security target tissue generate enough thermal damage and surrounding tissues almost not affected; A moderate amount of energy density security in the suitable time to provide strong enough energy output enough to damage the target tissue and normal tissue is almost not affected; Suitable skin protection measures to guarantee the target tissue enough damage, and skin almost not affected, so as to ensure the safety of the treatment.

## **Feature and Advantage**

- 1)Frozen environment , High-speed unhairing reduce pain truly minimally invasive .
- 2)Center wavelength of 808nm through the skin effect hair papilla directly and removal efficiency is higher than other hair removal device.
- 3)Use the 1800W power with the 500W semiconductor laser power supply high frequency and high efficiency light.
- 4)The compressor condenser that handle constant temperature cooling without intermittent work, more suitable for professional medical establishments.
- 5)Using the sliding and function handles ,continuous light sliding operation removal without from the site of action.

## **Application**

Permanent Hair removal of all body part : beard, lips,arm ,leg ,armpit hair,bikini line and so on .Improve pigment ,dark color of the mole ,burgundy birthmarks and other symptoms.

## Asia One 808nm Diode Laser



### 1. Theory of operating,

Semiconductor laser hair removal machine is freezing through technological innovation, so that the treatment of head to negative temperature, unlike other laser burns the skin, to better protect the skin, hair follicle melanin by specific wavelength of laser absorption, and the heat exchange to destroy hair follicle, the laser energy can be shot down into the hair follicles within the hair follicle is destroyed energy decomposition, so that the hair follicles to stop growing, because the laser energy about 15% of the skin can penetrate the dermis layer, the epidermis and more with melanoma of the hair follicle deep through the heat damage to hair follicle degeneration and atrophy; achieve permanent hair removal results.

**2. Main Features:**

1. Remove area: Lip hair, beard, chest hair, armpit hair, back hair, arm hair, leg hair, and unwanted hair outside the bikini line.
2. Secure high-speed repair damaged skin.
3. Whitening Firming ultra-efficient, ultra-soft skin can improve
4. Effective treatment for dry skin and rough skin
5. Effective treatment for large pores, fine wrinkles and redness.
6. Effective treatment freckles, reduced skin elasticity and dull skin, etc.

**3. Specifications**

| Items                   | Specs   |
|-------------------------|---|
| Laser wavelength / Type | 808nm / Diode laser                                 |
| Laser Energy output     | 0 ~ 40J / cm2 adjustable                            |
| Laser pulse width       | 1~400ms adjustable                                  |
| Pulse output frequency  | 1Hz   |
| Spot size               | 10*20mm   |
| Display interface       | 15" TFT Touch Screen                                |
| Contact cooling         | semiconductor cooling + water cooling + air cooling |
| Power Source            | AC 220V± 10% 50Hz                                   |
| Dimensions              | 420*380*1150mm                                      |
| Weight                  | 70 kg   |

**IPL**

**1. Theory of operating**

The system emits a broad spectrum of intense pulsed light through a small, smooth hand-piece. The light is placed over the skin gently. Cut-off filters <three different filters for option> in the hand-piece for change the wavelength range allowing it to be optimized for different applications and skin types. The light penetrates the tissue and is absorbed by either the blood when treating vascular lesions or the melanin when treating pigmented lesions and hair removal, and then damage them. The body's natural processes will remove the injured tissue, giving the skin a more even and youthful appearance. And at the same time, reproduction of the collagen in dermis will increase the flexibility of the skin and get rid of fine wrinkles.

**2. The function of IPL series**

1. Remove facial vascular lesions <Telangiectasis>
2. Clean and tighten visible pores,
3. Smooth fine wrinkles,
4. Make skin more flexible,

5. Improve the tone of skin and remove facial blemish,
6. Remove acne,
7. Removenunwanted hairs,
8. Skin whitening.

### 3. Specifications

|                     |                        |
|---------------------|------------------------|
| <b>Model</b>        |                        |
| <b>Wavelength</b>   | 430/530/640nm          |
| <b>Power</b>        | 1100w                  |
| <b>Pulse energy</b> | 1-50mj                 |
| <b>light spot</b>   | 10*30mm /<br>15*35mm   |
| <b>Dimension</b>    | 52*45*120cm            |
| <b>Net weight</b>   | 40 KG                  |
| <b>Voltage</b>      | 220V 50HZ<br>110V 60HZ |



#### Company advantages:

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

# IPL Device

## SH6.0 IPL Permanent Hair Removal Beauty Equipment



### Function

- Hair removal
- Spots removal: fleck, Sunburn, age pigment,chloasma etc
- Eliminate the wrinkle, perfect and smooth
- Breast enhancing
- Spider vein dispelled. Bottle nose, red spot treatment

- Shrink sweat pores, Ance, comido and acne ascars
- Activate the increase of collagen albumen and renew the elasticity of skin
- Improve the coarse skin on the body obviously.

#### **Specification**

- Light source: Intensive pulse light
- IPL energy: 10-50J/cm<sup>2</sup>
- IPL frequency: 1s,2s,3s,4s
- Pulse space: 5-60ms
- Pulse duration: 2-15ms
- IPL power source: 1600W
- Light spot area: 12\*30mm<sup>2</sup>/15\*35mm<sup>2</sup>/15\*50mm<sup>2</sup>
- Cooling system: semiconductor+ water+ wind+ cooling gel
- Input power : AC 220V(110V)/ 15A 50HZ(60HZ)
- Package Size: 53\*60\*133cm

#### **Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The engineers are the pioneers of optoelectronics in China, most professional and greatest innovation ability.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control

#### **Training:**

Once you buy machines from us,we supply professional training for you. Detailed user manual and online service are the basis training. Moreover, we teach how to use the machine and how to do the treatments step by step. Training certification will be sent to you after finishing training.

#### **After sale service:**

1 year warranty

Within 1 year warranty, we offer free repairing service and provide replacement parts for free.

1. we are the only one and the first one to supply 1 year warranty
2. we are very confidence about the quality of our products.
3. we try best to make clients work with the machine without any troubles after buying.

## SH6.1 Portable IPL Beauty Equipment



### Function

- Hair removal
- Spots removal: fleck, Sunburn, age pigment, chloasma etc
- Eliminate the wrinkle, perfect and smooth
- Breast enhancing
- Spider vein dispelled. Bottle nose, red spot treatment
- Shrink sweat pores, Ance, comido and acne ascars
- Activate the increase of collagen albumen and renew the elasticity of skin
- Improve the coarse skin on the body obviously.

### Specification

- Light source: Intensive pulse light
- IPL energy: 10-50J/cm<sup>2</sup>
- IPL frequency: 1s,2s,3s,4s



- Pulse space: 5-60ms
- Pulse duration: 2-15ms
- IPL power source: 1600W
- Light spot area: 12\*30mm<sup>2</sup>/15\*35mm<sup>2</sup>/15\*50mm<sup>2</sup>
- Cooling system: semiconductor+ water+ wind+ cooling gel
- Input power : AC 220V(110V)/ 15A 50HZ(60HZ)
- Package Size: 53\*60\*76cm

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The engineers are the pioneers of optoelectronics in China, most professional and greatest innovation ability.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.

**Training:**

Once you buy machines from us, we supply professional training for you. Detailed user manual and online service are the basis training. Moreover, we teach how to use the machine and how to do the treatments step by step. Training certification will be sent to you after finishing training.

**After sale service:**

1 year warranty

Within 1 year warranty, we offer free repairing service and provide replacement parts for free.

1. we are the only one and the first one to supply 1 year warranty
2. we are very confidence about the quality of our products.
3. we try best to make clients work with the machine without any troubles after buying.

## SH6.2 Portable IPL Photon Machine



### Function

- Hair removal
- Spots removal: fleck, Sunburn, age pigment, chloasma etc.
- Eliminate the wrinkle, perfect and smooth
- Breast enhancing
- Spider vein dispelled. Bottle nose, red spot treatment
- Shrink sweat pores, Ance, comido and acne ascars
- Activate the increase of collagen albumen and renew the elasticity of skin
- Improve the coarse skin on the body obviously.

### Specification

- Light source: Intensive pulse light
- IPL energy: 10-50J/cm<sup>2</sup>
- IPL frequency: 1s, 2s, 3s, 4s
- Pulse space: 5-60ms

- Pulse duration: 2-15ms
- IPL power source: 1600W
- Light spot area: 12\*30mm<sup>2</sup>/15\*35mm<sup>2</sup>/15\*50mm<sup>2</sup>
- Cooling system: semiconductor+ water+ wind+ cooling gel
- Input power : AC 220V(110V)/ 15A 50HZ(60HZ)
- Package Size: 53\*60\*76cm

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The engineers are the pioneers of optoelectronics in China, most professional and greatest innovation ability.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.

**Training:**

Once you buy machines from us, we supply professional training for you. Detailed user manual and online service are the basis training. Moreover, we teach how to use the machine and how to do the treatments step by step. Training certification will be sent to you after finishing training.

**After sale service:**

1 year warranty

Within 1 year warranty, we offer free repairing service and provide replacement parts for free.

1. we are the only one and the first one to supply 1 year warranty
2. we are very confidence about the quality of our products.
3. we try best to make clients work with the machine without any troubles after buying.

## SH888 Home Use IPL Device



### IPL Hair Removal Machine(iLight)

Home IPL hair removal machine is applied with intense pulse light technology. It brings a comfortable, safe, and permanent hair removal result. You can remove unwanted hair at your privacy of home.

#### Advantages of Deess IPL hair removal machine

1. Much faster and easier way of hair removal than other existing methods
2. Almost painless and less time-consuming compared to other treatments
3. Treats small areas of hair at a time

4. Less number of times for treatment
5. obvious and long lasting result etc
6. Permanent hair removal and hair reduction effect
7. You can feel the soft skin resulted from pore constriction after being treated.

**Where can be used?**

Lip hair, armpit hair, body hair on arms and legs, hair on places that affecting looks such as hair line

on forehead and bikini area.

**How long between two treatments?**

In the first 3 -5 treatments, do one treatment every 2 -3 weeks, in the following 6 -8 treatments, do one treatment every 4-6weeks.In the maintenance period, do one treatment every 2-3 monts for areas with hair regrowth.

**Featureas of the ipl hair removal machine**

- 1.Rechargeable,handheld,portable,convenient
- 2.Long lifetime:50,000 shots, each recharge can shot 200 times.

**Effects of the IPL hair removal machine**

Clinically proven, 90% hair is reduced after 7-10 treatments, a satisfied result can be gained generally . When being treated, most people experience a sense of warming.

**Treatment**

**Accessories: Glasses, warranty card, user manual and vedio.**

**How to use the ipl hair removal machine?**

**FAQ**

1. Can the IPL hair removal machine be used on the face, head or neck?

Recommend for non-facial sites treatment and below the neck:

on the legs, underarms, bikini line, back, chest, stomach, arms, hands, or feet. It has not yet been cleared in the U.S. for use on the face, head or neck.

2. Does the IPL hair removal machine really work?

Absolutely. IPL hair removal home use device is designed to gently disable hair growth so that your skin remains smooth and hair-free, for good.

3. When will I begin seeing results?

You will see immediately noticeable results, in addition, you will begin to see results after your third treatment and be virtually hair-free after ten. Be patient — the results are worth the wait!

4. How can I accelerate the results?

You will apparently see results faster if you have treatments twice a month for the first three months.

Following that, you still have to treat once a month for another four to five months to completely remove the hair.

5. Can the IPL hair removal machine be used on men?

Of course! We have already received many great cases, as men want the same permanent hair reduction as women.

6. Does it hurt?

Precisely speaking, the sensation varies by individual, but most people think the feeling as a light to medium rubber band snap on the skin, anyway, that feeling is significantly more comfortable than waxing,

Remember it is important to always use low energy settings for initial treatments.

7. Do I need to prepare my skin before using the IPL hair removal machine?

Yes. Start with a close shave and clean skin that's free of lotion, powder, and other treatment products.

8. Will hair grow back?

Yes, some of it will. However it will grow back in looking thinner and finer.

If you stop using IPL hair removal home use device hair growth may eventually return to its previous pattern.

9. Can I use it everyday?

It's not recommended to use everyday. Hair regrowth will not be sufficient for a successful treatment (1mm minimum length),. It is better to wait for at least 1mm of hair regrowth before doing next treatment.

10. Is the IPL hair removal machine likely to lead to a permanent hair removal?

As a matter of fact, there is no such thing as PERMANENT hair removal

(with the exception of electrolysis).

However, Clinical studies noted a marked reduction in hair growth

when IPL used for over 6 weeks.

11. Are there any side effects like bumps, pimples and redness?

Clinical studies show no lasting side effects associated with the proper use of

IPL hair removal home use device like bumps and pimples ;

However, people with hyper sensitive skin may experience temporary redness which fades within hours.

Applying smooth or cooling lotions after a treatment will help keep skin moisturized and healthy.

### **Warranty policy**

- All the beauty machines sold by us carry 1 year warranty and lifetime service. If you found your item to be defective, please return it to us in the original condition and package.
- Our warranty does not include replacement bulbs, any glass pieces, consumable parts such as filters or worn implements.

### **Payment term**

T/T, Western Union, Money Gram.

### **Delivery Time:**



2-3 days after payment

**Packing Size,**

22cm\*15cm\*8cm

**Gross Weight:**


2kgs





# Infrared Pressotherapy Machine

## SH5.0 Infrared Pressotherapy Slimming Machine

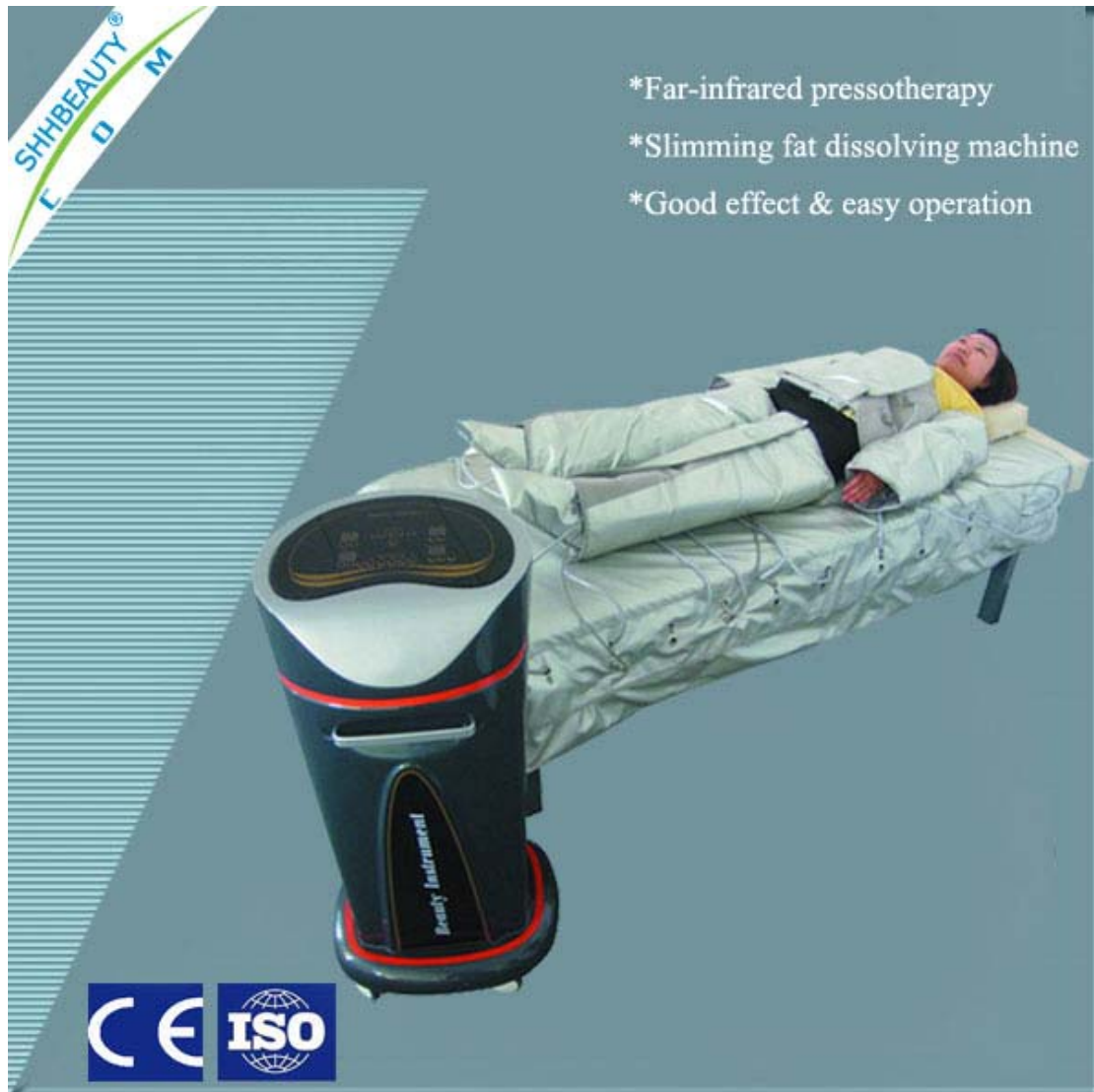

### 1. Technical specification

Material:ABS

Output power:350w

Warranty:3 years

Power supply:220v/110v, 50hz/60hz

Package size: 56\*61\*111cm

GW: 40kgs

2. Core Technology

Pressotherapy

Infrared

3. Function

body slimming


blood circulation promotion

4. Main accessories

1pcs body suit for body slimming and circulation of blood promotion

Company advantages:

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

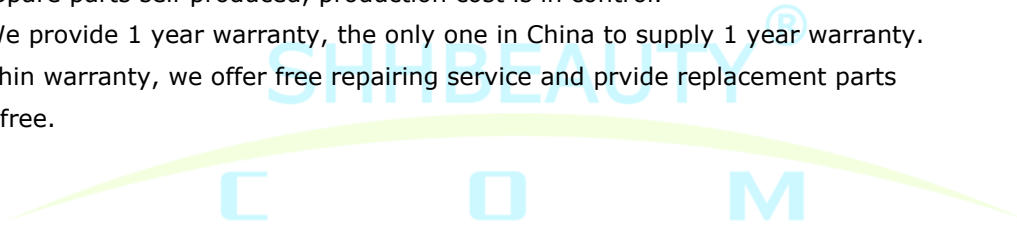

## SH5.1 Air Pressure Pressotherapy Infrared Equipment



### 1. Technical specification

Material:ABS

Output power:350w

Warranty:3 years

Power supply:220v/110v, 50hz/60hz

Package size: 55\*46\*103cm

GW: 40kgs

### 2. Core Technology

Pressotherapy

Infrared

### 3. Function

body slimming

blood circulation promotion

### 4. Main accessories


1pcs body suit for body slimming and circulation of blood promotion

### Company advantages:

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.



## SH5.2 Pressotherapy Air Pressure Slimming Equipment

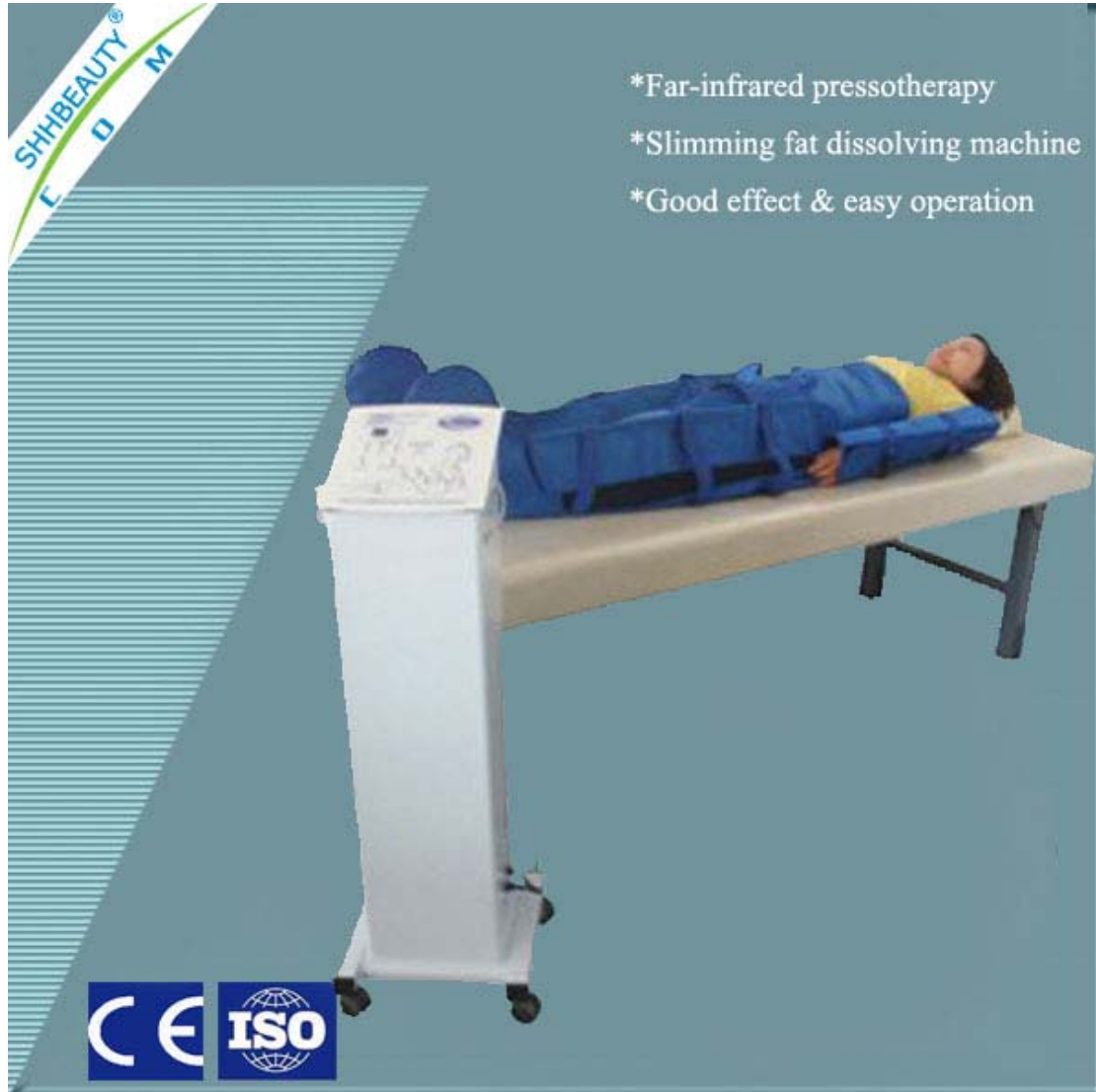

### Specification,

Power supply, 220V/110V, 50Hz/60Hz

Output power, 750W

Gross Weight, 30kgs

Packing, 47\*52\*110cm

### Features of Air pressure slimming cloth:

1. This equipment applies air pressure technology.
2. Unique design with zipper, size adjustable according to different figure.
3. easy to control, fit to household use.
4. Effective to promote blood circulation, lymph drainage, softens cellulite, eliminate fatigue, and relive pain.

**FUNCTIONS:**

Weightn loss

Slimming

Body building

**ACCESSORIES:**

Air-pressure machine

Slimming clothes


Both the clothes and the machine connected by the power lines attached on the outputs

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.



## SH5.4 Portable Air Pressure and Infrared Equipment

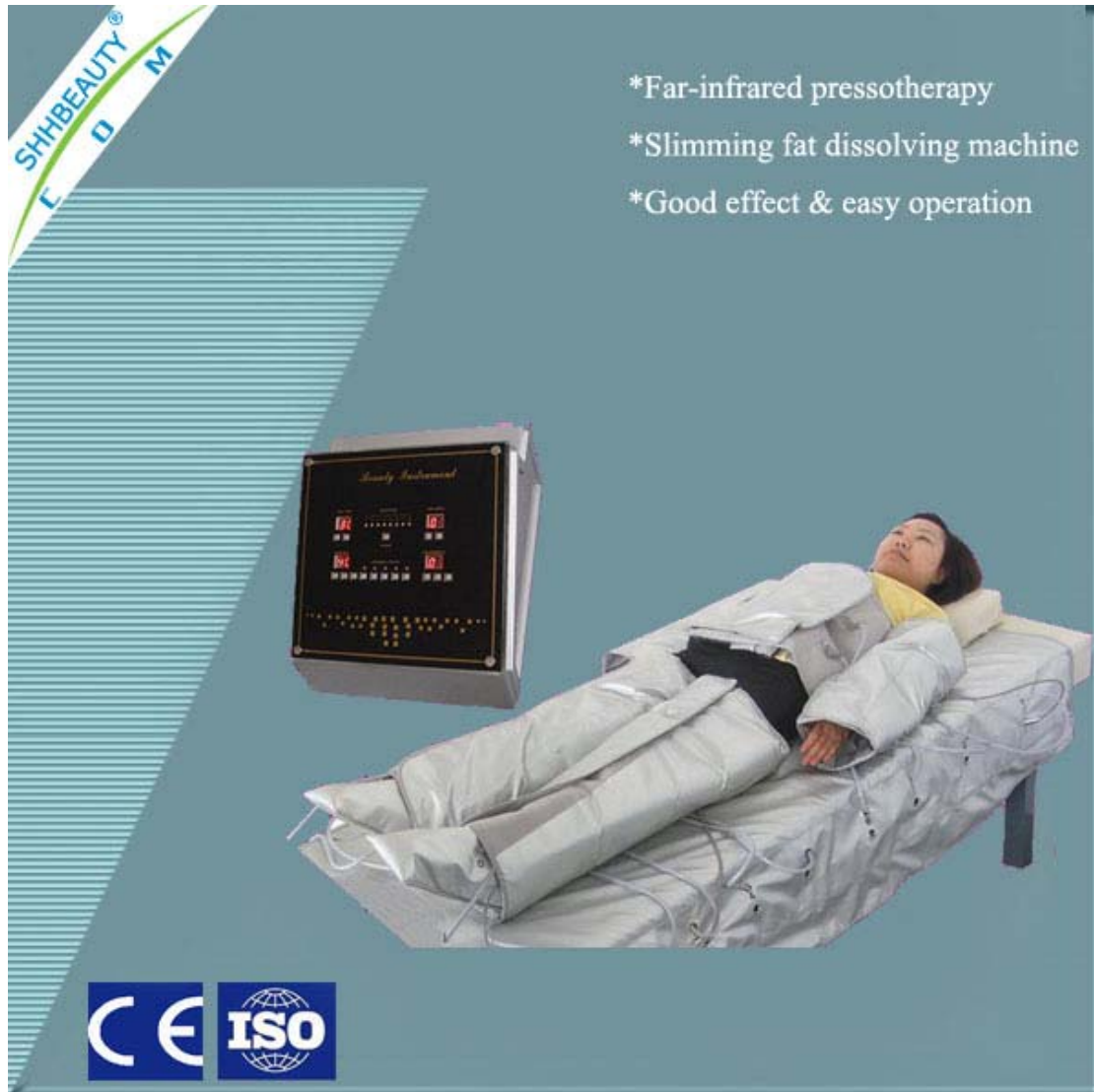

### SPECIFICATION,

|                                |             |
|--------------------------------|-------------|
| Package size                   | 97*55*42cm  |
| Power                          | 70W         |
| Color                          | white       |
| Voltage                        | 220V/110V   |
| W.G                            | 45kg        |
| Outputs                        | 4           |
| Control panel                  | Button-type |
| Blanket size                   | 170*77cm    |
| Intensity of pressure air wave | 80 KP       |
| Warranty                       | 1 year      |
| Temperature range              | 0-69°C      |

**Theory information:**

Air vibrant massager aims to promote normal flow of lymphoid liquid and sweep the body rubbish. Unlike vibrating by adhesive piece before, it works by means of pressure depend on envelopes. With assist of the apparatus, beauticians adopt several basic techniques as circling and steering massage, so it can improve blood circulation and toxin excreted to make capillary passable. Air De-toxin system ( lymph edema compression therapy system) promotes the circulation of blood and lymph and make recovered patients leg or arm having various edema or troubles by stimulating the muscular organization by repeating the contraction expansion program that pumps and deflates air through each chamber of the cuff that exerts the pressure from distal to proximal of body gradually. Air De-toxin system has many models that have the various function and the cuffs (pants, leg, arm, shoulder, skirt) having 2-chamber ,4-chamber, 6-chamber or 12-chamber. According to each model, air de-toxin system is a physical therapy appliance, which cures the patients with the unique program. FEATURES: As it massages with smooth air pressure, there is no danger from any electric shock or electromagnetic wave, and just give nice and fresh feeling. It is made pressure and time of massage to be controlled according to user constitution and symptoms so that whoever can use safely DR life It is very portable to be used freely and easily. It can be used partially (arm, leg, waist) for your useful anyway. It has the gauge that can check the pressure.

**Function of infrared blanket:**

1. Remove unwanted fat drainage toxin
2. Accelerate matabolism, strengthen nutrition using and dispel tire
3. Do physiotherapy on muscles and neuropathic diseases, eliminate liver strain and unwanted horde pains
4. Enhance immune function, prevent cancer cells dilation
5. Diminish pain and inflammation

**Features of Air pressure slimming cloth:**

1. This equipment applies air pressure technology.
2. Unique design with zipper, size adjustable according to different figure.
3. easy to control, fit to household use.
4. Effective to promote blood circulation, lymph drainage, softens cellulite, eliminate fatigue, and relive pain.

**CORE TECHNOLOGY:**

Air-pressure

Infrared slimming

**FUNCTIONS:**

Weight loss

Slimming

Body building




**ACCESSORIES:**

Air-pressure machine

Slimming clothes

Both the clothes and the machine connected by the power lines attached on the outputs.

## SH860 Infrared Pressotherapy Slimming Machine

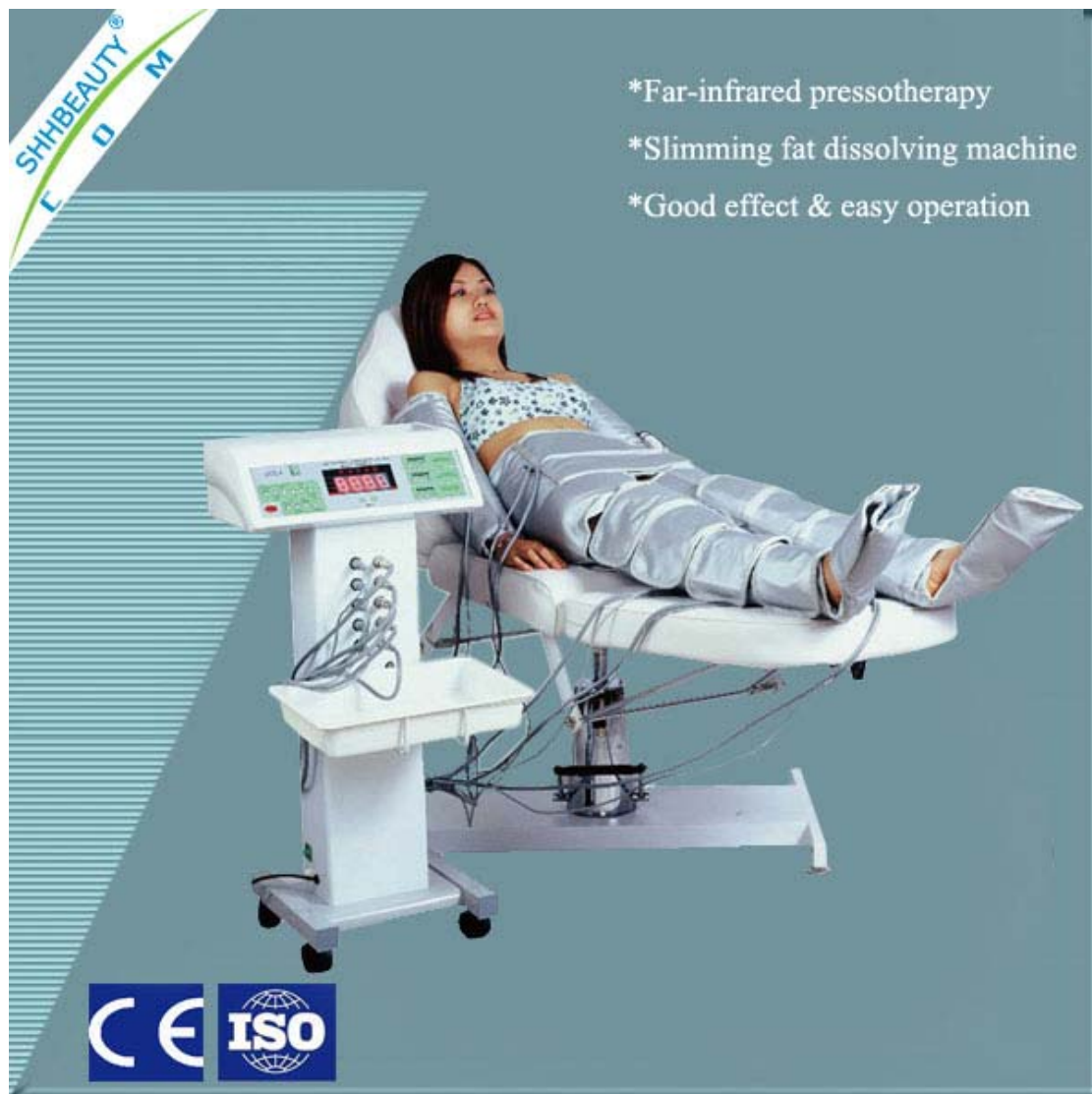

**Theory:**

Perfectly combining the functions of removing lymph & ejecting fatness, dissolving fatness, breaking fatness, magnetic therapy and fasten skin, also gathering the exercises of slimming, body-beauty and massaging.

Emitting 32 signal with different wavelength and inputting into the body to produce 32 moving methods and dissolve the fatness, meanwhile, it combines Japan fibre and emit 8UM-14UM ultra red dissolve the fatness in the deep body, exhaust the fatness through lymph strong movement, with its specially-designed strong magnetic and applying with

super-length wave, while the movement of the fatness, activate the cell and make the muscle bright and smooth, to achieve the result of slimming, body beauty and massage.

**Function:**

1. FAR-infrared dissolve fat.
2. Super-long-wave break fat.
3. Vent Lymphocyte and fat.
4. Magnetic Therapeutics and shrink skin all the way.

**Packing & Delivery:**

1. packing Size:96.8X59.4X47.5cm
2. G.W:26.4kg
3. N.W:21kg
  
4. Material:ABS
5. Output power:350w
6. Power supply:220v/110v, 50hz/60hz

**Warranty:** 12 month

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.

## RF Radio Frequency

### SH7.2 Latest Portable Wrinkle Removal Facial Lifting RF Beauty Machine for Home Use.



#### How does RF System Work ?

Bipolar RF technology offers new scientific development in cellulite treatment, body sculpting, anti-aging and skin tightening. It involves applying electromagnetic (non-ionizing) waves to the body in order to return physiological processes that have slowed down to normal. These waves carry heat and energy deep into the skin creating an endothermic effect that has a positive effect on the lymphatic and vascular systems and on the adipose-connective structures.

## Functions

Bipolar RF waveform with cold compress and four wave principles, Accelerate blood circulation, activate cells, and strengthen the metabolism, a fundamental solution to the skin problem.

- 1) Sculpture facial continuers, such as: facial slimming and chin lose weight, diminish wrinkles, laugh when the wrinkles, contractiorelaxation of the skin, the use of the effect to replace and exceed the effects of surgery lift.
- 2) Lymphatic dredge and remove excess water and adriamycin, eliminate the bags under the eyes, black eyes shrink pores, treat acne.
- 3) Soften skin and accelerate the decomposition of melanin, whitening dispel macular.
- 4) Systemic anti-aging, neck and back contraction modification, shaping around the waist and abdomen, shape the buttocks and legs.

## Specifications

Input voltage: AC100-240v

Frequency: 50/60HZ

Output voltage: DC9V-2000Ma

Warm temperature: 40-60centigrade/104-140 °F

Cold temperature: 0-7centigrade/32-44.6°F

Packing Size, 34\*24\*23cm

Packing included

- Main Unit
- Handles:
  - RF needles probe airwaves x 4
  - RF semicircles x 2
  - frozen machine x 1

## RF3.6 Multipolar RF Vaccum Dissolving Fat System



Vaccum BODY SLIMMING machine can promote tissue metabolism, repel the cellulite, tighten the skin, strengthen the skin elasticity and so on. Its effect will last for a long time. Meanwhile, it can enhance the elasticity of muscle and improve the metabolism of cells.

The Tripolar Radio-frequency (RF) utilizes 10MHz/S electric wave, which can permeate the epidermis to take effect on the underlying tissue. High-frequency shake will make resistance of the underlying tissue move and produce high energy, which makes the temperature of the underlying tissue rise instantly. Since high heat causes immediate dissolve of fat and the body become slim, hence the desirable effects of skin-tightening, wrinkle-removal and figure-shaping, etc. This machine working with cavitation and RF, it can help to remove our fat more effectively.

**Function**

Tripolar RF for eyes;  
Tripolar RF for face;  
Seven Polar RF for body;  
Eight polar RF with Vaccum for Body Slimming;  
Touch Screen

**Parameter**

Operating Voltage: AC220V/50Hz  
Power: ≤50W  
NW: 45KG,  
GW: 65KG  
RF FREQUENCY: 0.3M/1MHZ/2MHZ  
Vaccum pressure: 0.1MPA  
Packing: 77cm\*48cm\*125cm

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The engineers are the pioneers of optoelectronics in China, most professional and greatest innovation ability.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control

**After sale service:**

- 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.
1. we are the only one and the first one to supply 1 year warranty
  2. we are very confidence about the quality of our products.
  3. we try best to make clients work with the machine without any troubles after buying.

## RF395 RF Electric Wave Tighten Skin Instrument



### Theory

wrinkle removal system makes use of E light and RF technology, providing with ideal temperature scatter for deep dermis layer, helpful for the exact treatment for targeted tissue, effective for deep dermis and connective tissues. It will stimulate collagen and make it rebirth without affecting the formal skin. It will safely remove wrinkle and improve skin texture.

### Main effects:

- a. resume skin youth and skin tighten physical therapy
- b. acne physical therapy
- c. special eyes physiphysical therapy ,dark cycle
- d. lymph drainage physical therapy
- e. lip physical therapy

**Specification**

korean RF operating head  
monopole design  
foot controller  
451khz  
power < 200w

Packing Size: 46\*36\*31cm

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.





# Oxygen Injection Instrument

## G228A Omnipotent Oxygen Jet Therapy Machine



### **Detailed Description,**

The omnipotent skin oxygen-injection instrument applies the world's advanced PSA theory. It can separate medical-use purified oxygen (high density, no dust, no bacteria) from air in normal temperature; once the power supply is connected, it can absorb oxygen, inject oxygen and spurt oxygen just with one machine. Therefore, it can solve the difficulties of anoxia. Inject 98% purified oxygen and active nutrition distillate to the deep-seated of skin through impulse force, in order to accelerate the metabolism of cells, whiten and tender the skin.

**Specification,**

Inject oxygen  
 Spurt oxygen  
 Inhale oxygen  
 Photon Skin activation  
 BIO crinkles elimination  
 Lymph smooth  
 MEAS.56cm\*49cm\*136cm  
 42.5cm\*39.5cm\*123cm  
 GW:52KG

**Applications:**

1. Skin rejuvenation, remove facial blemish, tighten large hair pores, deep cleaning, and improve the flexibility and tone of skin.
2. Scar removal; Remove scar caused by laser treatment, burn, surgical etc.
3. Wrinkle removal; Remove fine lines, wrinkles and spider veins.
4. Acne removal; Remove acne, comedones, blackheads and acne scar.
5. Restrain the inflammation.

**Feature and benefits:**

1. New idea for skin rejuvenation, safe treatment with good result.
2. A vast range of applications; Many in one machine.
3. Safe and effective for all skin types.
4. Friendly design, easy to operate.
5. Skin scan system involved, before after picture can be taken and displayed.
6. Air compressor inside, reserved oxygen bottle connection port makes two options available.
7. emitting compressed air (using air compressor inside); Or emitting medical oxygen (connected with oxygen bottle when operating).
8. Skin Scanner and LCD Screen involved make easy contrast.

**Therapeutic heads,**

1. Expelling of toxin, it can help you expel the toxin of body;
2. Two inject oxygen therapeutic heads;
3. One Spurt oxygen therapeutics, it is can whiten and tender skin etc.
4. Soft photon therapeutic head with the function activating skin, accelerating the blood circulation, promote the metabolism etc.
5. BIO crinkles elimination rods. It can be penetrate into dermis directly and release the energy, which can penetrate into skin from inner to out layer by layer, It has an effective solution to the problem of skin such as pachylosis and speckle caused by vitamin uneven or climate change, so that it can make skin bright, delicate or slender, shiny.

## G882A Hyperbaric Oxygen System



### Specification:

Input: AC220v/50hz AC110v/60hz

Power: 450W

Oxygen purity: 96.4%

RF frequency: 5MHZ 75KHZ

Wavelength: 330-990nm

Music Output: Self-selection

Man- machine control: Color Multimedia and large touch screen

Packing Size: 53\*59\*148cm

Gross Weight: 52kgs

## **HYPERBARIC OXYGEN SYSTEM HBO**

1. HYPERBARIC LIGHT MASK, oxygen, LED, music therapy
2. MULTI POLAR RADIO FREQUENCY , bipolar rf ,tripolar rf
3. COMBINE COLD LIGHT, MUSICOTHERAPY,RF, OXYGEN INJECT, OXYGEN SPARY,
- 4 magic hands: micro current therapy,

## **BIO AND SUPERSONIC FUNCTIONS**

### **FUNCTION:**

SKIN OXYGEN THERAPY, STIMULATE SKIN REGENERATION, BLOOD CIRCULACION  
IMPROVE BODY METABOLISM, AND COLAGEN REPRODUCTION,

### **ADVANTAGES:**

- 1 ELECTRONIC OXYGEN GENERATION
- 2 PATENT HYPERBARIC MASK WITH MUSICOTHERAPY
- 3 MULTI FUNCTION

### **HANDLES:**

- 1.OXYGEN SPARY GUN
2. OXYGEN INJECT HANDLE
3. OXYGEN MASK WITH LIGHT, AND MUSICOTHERAPY
4. SUPERSONIC HANDLE ,2 PIECES
- 5 BIPOLAR RF FOR EYE
- 6 TRIPOLAR RF FOR FACE

### **How is about the design of this machine?**

The design and the current is fully comply with European electrical safety standards, absolutely safe to human body, the built-in anti-electromagnetic interference devices, can effectively protect against electrical and magnetic interference, maintain the current stability, to ensure that procedures effective and stable.

### **Are you real manufacturer or just trade company?**

We are one of the leading real beauty equipment manufacturer with more than 200 employees. We have more than 15 years experience for production of beauty machines. We have our own strong designing and developing new products ability in beauty industry.

### **Could you introduce your after-sales service?**

We provide customers 1 year warranty time and life-time technical support. Within the warranty time, customers could get free accessories for repairing the machines purchased. We have a professional after-sales team consists of 10 beauty machine engineers. Also our beauticians are so glad to help customers for anything to get good effectiveness of beauty machines treatment.

## PDT LED Soft Photon Machine

### SK8 PDT LED Soft Photon Light Therapy Instrument



**Specification:**

Voltage: AC110V/AV220V

Light: Red/Yellow/Green/Blue/Purple/Radiation

Wave length: 200NM-900NM

Out frequency: 1-30Hz

GW: 24KG

MEAS. Of package: 76cm\*66cm\*116cm

**Function:**

Improve symptoms such as aging and loosening skin, coarse pores, slender wrinkles.

- 1,improve pigmentary pathological changes,such as freckle, sunburn,senile plaques;
- 2,improve dark complecion caused by bad metabolism or poor circulation.
- 3,repair and nurse the damaged skin.
- 4,do detumescence,inflammation diminishing and scar elimination treatment for oil acne effectively.

**Theory,**

Photodynamic is cold light, and could not produce high temperature, and could not burnt skin. Transfer light energy into cell internal energy, accelerate the process of cell growth cycle, stimulate the production of collagen fiber cell protein. It can cause unexpected effects of repair aging skin, acne skin, skin patches, wrinkled skin, allergies and the sun burns the skin, what's more, which is could not bring any harm or uncomfortable to the skin.PDT is the world's top scientific and technology achievement, which applies the response of photodynamic to skin treatment and maintenance, It is a facial skin care expert, with simple operation and no side effect.

# Fractional co2 Laser

## SH9.0 Fractional Co2 Laser Equipment



### Technical parameter

- 1.Laser type:10600nm,CO2 laser
- 2.Optical maser:RF laser
- 3.Operating system :Korean auto intelligent controller system
- 4.Output pattern:Circle, square,ellipse , sexangle, rectangle, beeline, etc. for your option
- 5.Spot diameter:0.12-1.25mm
- 6.Laser power:10W/30W
- 7.Display screen:8.4 inch colorful touch screen
8. Light transfer system:7 Joint hinged arm
- 9.Guide light source:Semiconductor red light

10. Dot interval:0.5, 1.0mm
- 11.Treatment depth:0.1-2mm
- 12.Working volt:220V/50HZ 110V/60HZ
13. Dimension:53cm\*38cm\*139cm

### **Theory of operating**

Ultra pulse treatment is a "fractional" treatment that produces thousands of deep, tiny columns in your skin, called microscopic treatment zones. FLT is "fractional" in that it intensively treats old epidermal pigment cells without affecting any of the surrounding tissue, which reduces healing time about 24 hours) and promotes the skin's natural process for creating new, healthier tissue to replace the damaged cells. Ultra pulse treatment is a safe, non-invasive procedure that allows you to rediscover the fresh, healthy skin of your youth.

### **The scope of treatment**

- 1.Acne scars
- 2.Actinic keratosis
- 3.Basal cell carcinoma
- 4.Burn debridement
- 5.Laser skin resurfacing
- 6.Eliminate the scar Traumatismoperation
- 7.Improve texture tone and pore size
- 8.Erase unwanted brown spots
- 9.Smooth wrinkles around the eyes and mouth
- 10.Eliminate striae gravidarum
- 11.Rejuvenate your neck, chest and handsets



# Breast Enhancer Mamchine

## M7 Super Breast Enhancer Machine



### Detailed Product Description

1. breast enhancer
2. Electro-Stimulationsystem
3. LED Photomodulation Mask
4. COSCO IRnipple rejuvenation system
5. CE

**Feature**

1. Electro-Stimulation system
2. LED Photomodulation Mask
3. COSCO IRnipple rejuvenation system
4. Infrared heat activation system
5. Vacuum suction fat removal system
6. Rhythmic aerobic vacuum system
7. Synchronous rhythm of heat rate
8. Radio frequency thermal system
9. Big color touch-screen with multimedia, control and display treatment procedures.
10. USB multimedia data transmission
11. The design and the current is fully comply with European electrical safety standards, absolutely safe to human body, the built-in anti-electromagnetic interference devices, can effectively protect against electrical and magnetic interference, maintain the current stability, to ensure that procedures effective and stable.

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.

## N1 Vacuum Breast Enlargement



### Function:

1. Breast enhancement,
2. Nipple care,
3. Mold figure, body shaping
4. Vacuum
5. Lymph stimulation

### Features:

Sucking

Make the operation head keep good contact with skin and pinch it up quickly

Rhythm Massage

Make the operation head keep good contact with skin and do Rhythm massage as pre-programmed

Winding Massage

Make the operation head keep good contact with skin and slide windingly as shown in the picture

Sliding

Make the operation head keep good contact with skin and slide smoothly to relax.

Remove Fat

Make the operation head keep good contact with skin. Slide as shown in picture and pinch the skin quickly to remove fat

Parameter

**Dimension:** 66×37×26(CM)

**G.W.:** 9.5KG

**Power:** 45W

Frequency: 50/60Hz

Voltage: 110V/220V

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, the only one in China to supply 1 year warranty. Within warranty, we offer free repairing service and provide replacement parts for free.



# Micro Needle Roller

## BIO50 Micro Needle Roller



### Detailed Product Description

Micro-needle for skin care  
 stainlesspin and titanium alloy needle  
 five kinds of needle lenght,  
 two kinds of package

### Theory of Microneedle Therapy:

Microneedle therapy uses a knead rod inlaid with 192 slender steel needles, accompanying with functional nutrition liquid, regularly rolling in the places which need to cure, it can

puncture more than 300, 000 micro vessels on the epidermal within five minutes, form a efficient nutrient delivery system in the epidermis and subcutaneous tissue, with which the active ingredient of the nutrition can entry into skin effectively. Then make skin maintenance, skin care, skin treatment and skin beauty to be perfect. Microneedle therapy has the advantages of safe, reliable and multifunction. It can be used to beauty; Body slimming; Hair regeneration, etc. It comes to be the preferred method for skin anti-aging.

Function of the product: Stimulate skin to absorb the medicine and other active ingredient efficiently. Needle material is stainless steel. Characteristic of Micro-needles therapy

Use natural healing ability of human body, not stiffen skin or make skin cells shrinked.

Microneedles can make a lot tiny channels and thus make ingredient effectively penetrate the skin, which can produce more than 300,000 small holes in just 6 minutes, so as to making direct effects in the dermis.

Minimal side effects.

Micro-needles therapy is effective, safe and reliable with simple operation and reasonable price.

The Microneedle care is simple, can do common outdoor activities, including basking, but please note anti-sunburn supplies are needed.

**Applied range:**

1. Scar removal including acne scar removal or treatment.
2. Stretch mark removal
3. Anti ageing.
4. Anti wrinkle
5. Cellulite treatment/ cellulite reduction or removal.
6. Hair loss treatment/ hair restoration
7. Hyper pigmentation treatment.

**Spec.: 0.2mm, 0.5mm, 1.0mm, 1.5mm, 2.0mm**

0.2 mm

Skin Restoration,  
Improve Alopecia,  
Homecare Therapy,  
Active ingredient delivery

0.5mm

Skin Whitening,  
Improve skin pigmentation,  
Improve wrinkle (eye area),  
Improve Alopecia,  
Collagen Induction,  
Drug delivery

1.0 mm

Improve Wrinkle,  
Collagen Induction,  
Improve Skin Scar,  
Cellulite treatment,

Skin Thickness,  
Face Firming,  
Face Body Lifting

1.5 mm  
Pregnancy lines  
Thick wrinkle  
Large pore

2.0mm  
Dark sore concavex hole  
Pimple scar  
Severe decline of elasticity

## Microdermabrasion Machine

### NV-60 Diamond Jet Peel Machine



**Functions:**

Diamond Dermabrasion provided a non-surgical skin refinish procedure, by using sterile diamond heads to abrade or rub off the top skin layer, then vacuuming out the particles along with any dirt and dead skin back up. This procedure removes skin debris, imperfections, blemishes, wrinkles and unwanted pigmentation on the skin. Moderate usage of Diamond Dermabrasion, accompany with skin products that penetrate through top skin layer that reach the dermis level, helps to replenish natural nutrients, restore cell activity, and promote skin health. The results after Micro Dermabrasion treatment is usually a healthy, glowing, beautiful skin.

**Specifications:**

Packing Size: L37.5\* W36.5\* H29.5(cm)

G.W. 5KGS

Power: 45W

Voltage: 110V/ 220V

**Package Included:**

- 1) Diamond Dermabrasion Machine X 1
- 2) PVC vacuum hose 6\*4 X 1
- 3) Stainless steel wands X 2
- 4) Diamond tips X 9
- 5) Cotton filter X 1 package
- 6) Additional connector, tube rings and fuses
- 7) Full detailed English manual/instructions

**Package:**

Standard Export Carton with fully Polyfoam or Pearl Wool

**Company advantages:**

1. Specialized in production of medical and aesthetic equipments for many years.
2. The professional and innovative engineers are the pioneers of optoelectronics in China.
3. ISO9001 proved Quality control and test system, each machine is in the best quality control.
4. CE proved machines, the key for you to legally use and sell the machine.
5. Spare parts self produced, production cost is in control.
6. We provide 1 year warranty, Within warranty, we offer free repairing service and provide replacement parts for free.



## NV1005 Micro Current Facelift Machine



### Specification:

|              |                               |
|--------------|-------------------------------|
| Packing size | 58*49*38cm                    |
| Power        | 13.5W                         |
| Gross Weight | 4.6kg                         |
| Voltage      | 220V-240V/50Hz,100V-120V/60Hz |

### Can be used on:

Micro wave Anti-wrinkle massage suitable for sagging, ageing, excrescence skin and eye bag muscle.

Skin Maintain the pores and vessels. Clean up pores and vessles. Close up pores tingten nuckles. Improve nutrition absorption and permeability. at the same time, prevent aging and help face reforming.

Body Purify our body and blood. Smash fat cells and burn the fat have a same effect of arranging then makes our body look more airy-fairy.

**Struction:**

- 1.Turn on the power.
- 2.Insent the two handle.
- 3.Select the time as you like about 15 minutes or 30 minutes.
- 4.Select the strong or weak intensity.
- 5.Select the mode just you like sach as hip,bust,neck or face.

**Functions:**

- 1.Tone and condition facial muscle
- 2.Smooth laughter lines, wrinkles and fine lines
- 3.Tighten skin and improve elasticity
- 4.Boost circulation and increase supplies of nutrients to skin
- 5.Body sculpting and toning
- 6.Lymphatic drainage
- 7.Minimize eye bags and dark circles
- 8.Improve skin complexion

**After sale service:**

1 year warranty, within warranty, we offer free repairing service and provide replacement parts for free.



## Contact US

Shenghua Beauty

Tel: 0086-310-7023715

Mobile, 0086-15830761968

Fax: 0086-310-7023715

E-mail: [salesmanager@shhbeauty.com](mailto:salesmanager@shhbeauty.com), [shenghua@shhbeauty.com](mailto:shenghua@shhbeauty.com)

MSN: [shhbeauty@hotmail.com](mailto:shhbeauty@hotmail.com)

Website: <http://www.shhbeauty.com>

Address: No.485, Zhonghua North Street, Handan City 056004, Hebei Province, China